

Manual Wheelchairs

C O L L E C T I O N



karma

A Better Wheelchair, A Better Fit



Karma


"I feel alive only when wheels are spinning..."

These words came from a wheelchair user express the importance of mobility and that's the importance we are committed to support. Because these brave people deserve better wheelchair that suits their needs.

Established in 1987, "Fit the Needs of the Users" is always the core value that motivate Karma Medical keep moving forwards. The three major concepts of our innovation are: product development with specific features to correspond to users' medical needs, cutting edge ergonomic design and environmentally-friendly. We provide the best wheelchairs and help each customer start a new life with freedom and comfort.

In the past 30 years, we have built up over 4,000 dealer services channel in more than 40 countries successfully, supporting millions of wheelchair users and their families, enriching the lives of physically challenged people and sharing joyful and touching moments all around the world.

Global Marketing and Service Network

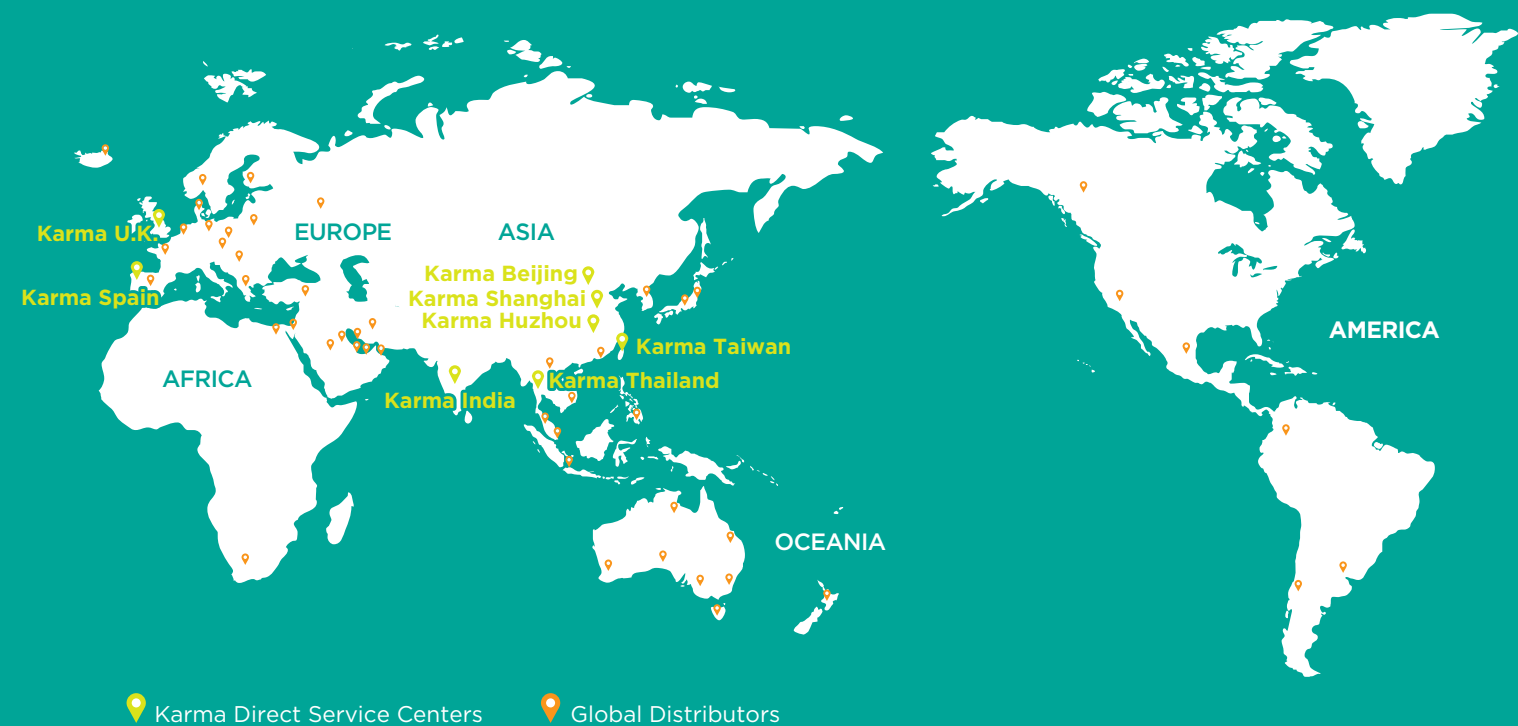

EUROPE

- U.K.
- Spain
- Belarus
- Benelux
- Czech
- Denmark
- Finland
- France
- Germany
- Greece
- Iceland
- Norway
- Poland
- Russia
- Spain
- Slovakia

ASIA

- Beijing
- Shanghai
- Huzhou
- Taiwan
- Thailand
- India
- Hong Kong
- Indonesia
- Iran
- Israel
- Japan
- Korea
- Malaysia
- Philippines
- Singapore
- Thailand
- Vietnam

AMERICA

- Argentina
- Canada
- Colombia
- Chile
- Mexico
- U.S.

OCEANIA


- Australia
- New Zealand

WEST ASIA

- Bahrain
- Kuwait
- Oman
- Qatar
- Saudi Arabia
- Turkey
- UAE

AFRICA

- Egypt
- South Africa

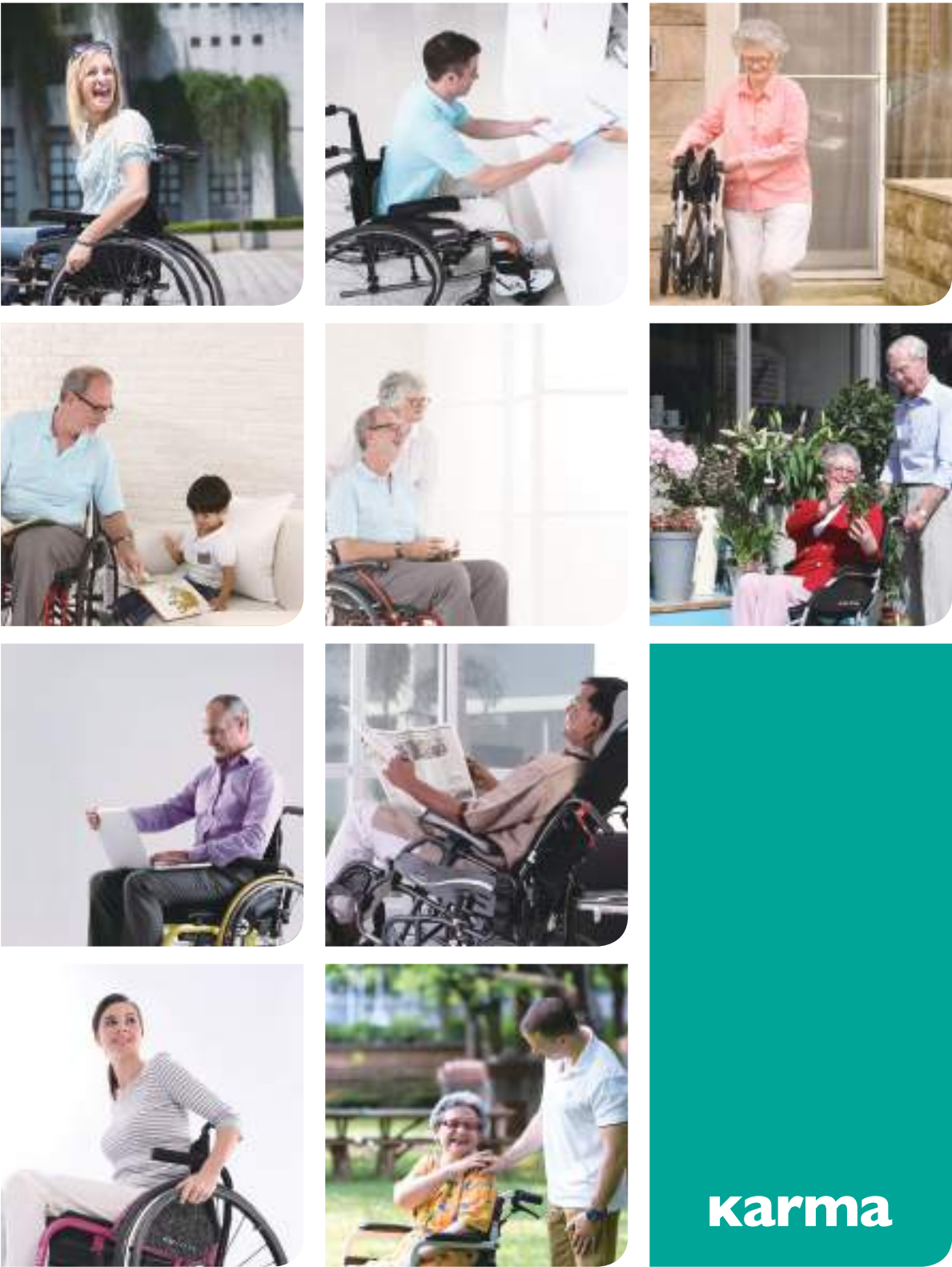

karma

Selecting a right wheelchair is crucial to overall quality of life. A wheelchair must be sized properly and include the appropriate features to meet personal lifestyle, body condition and medical needs. Please follow the simple steps in the guide below to have more correct concepts about how to choose a wheelchair.

1. Fit your Lifestyle


Where and how will the wheelchair be used? For a frequent user, a wheelchair with detachable or flip-back armrests will suit your transfer requirement. For a traveler, a lightweight, foldable and compact wheelchair is suggested. A compact wheelchair with excellent maneuverability is recommended in limited space. Furthermore, a suitable rear wheel size needs to consider the ability to propel.



2. Fit your Body Dimensions

a. Hipbone area

As a high risk group of pressure ulcers, a daily user should choose the wheelchair with ergonomic seat to relieve the pressure effectively. For occasional or temporary use, a basic wheelchair could be sufficient, yet the user is able to reposition body from time to time to relieve pressure.

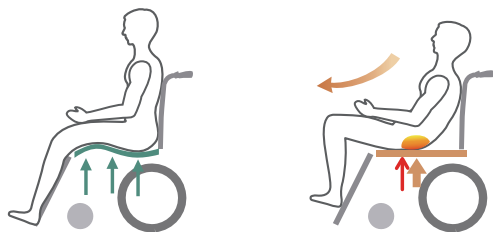


Ergonomic Seat Wheelchair

Standard Wheelchair

b. Back

“Reclining function” is ideal for the user with limited mobility to sit in the wheelchair comfortably. Also, it helps the attendant take good care of the user more easily. A wheelchair with “tilt-in-space function” is able to provide significant pressure relief and ensure good postural support.

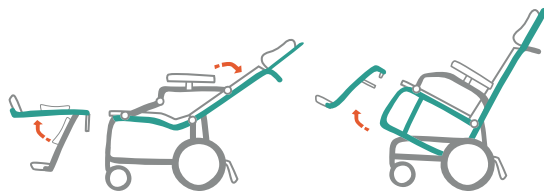

Recline Wheelchair

Tilt-in-space Wheelchair

c. Seat width

When a user sits in a wheelchair, the space from legs to side panels should be at least two fingers wide. This will ensure the user is seated comfortably and adequately supported.




d. Seat depth

How to make sure a proper seat depth? Subtract 2 inches from the length of backside of the knee to lower back. This will ensure the user is well supported and the seat upholstery will not rub against legs.



e. Footplates height

The feet should be placed properly on the footplates and the knees should be in line with the hip.

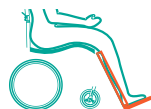


3. Fit your Medical Needs

The medical needs have to be considered when you choose the wheelchair. You can consult your own therapist to have more suggestions.

What is Karma FIT ?

We believe every wheelchair is not identical. Wheelchair users should be able to find one that fits individual needs. The wheelchair does not merely meet their basic requirements but actually improve their quality of life.

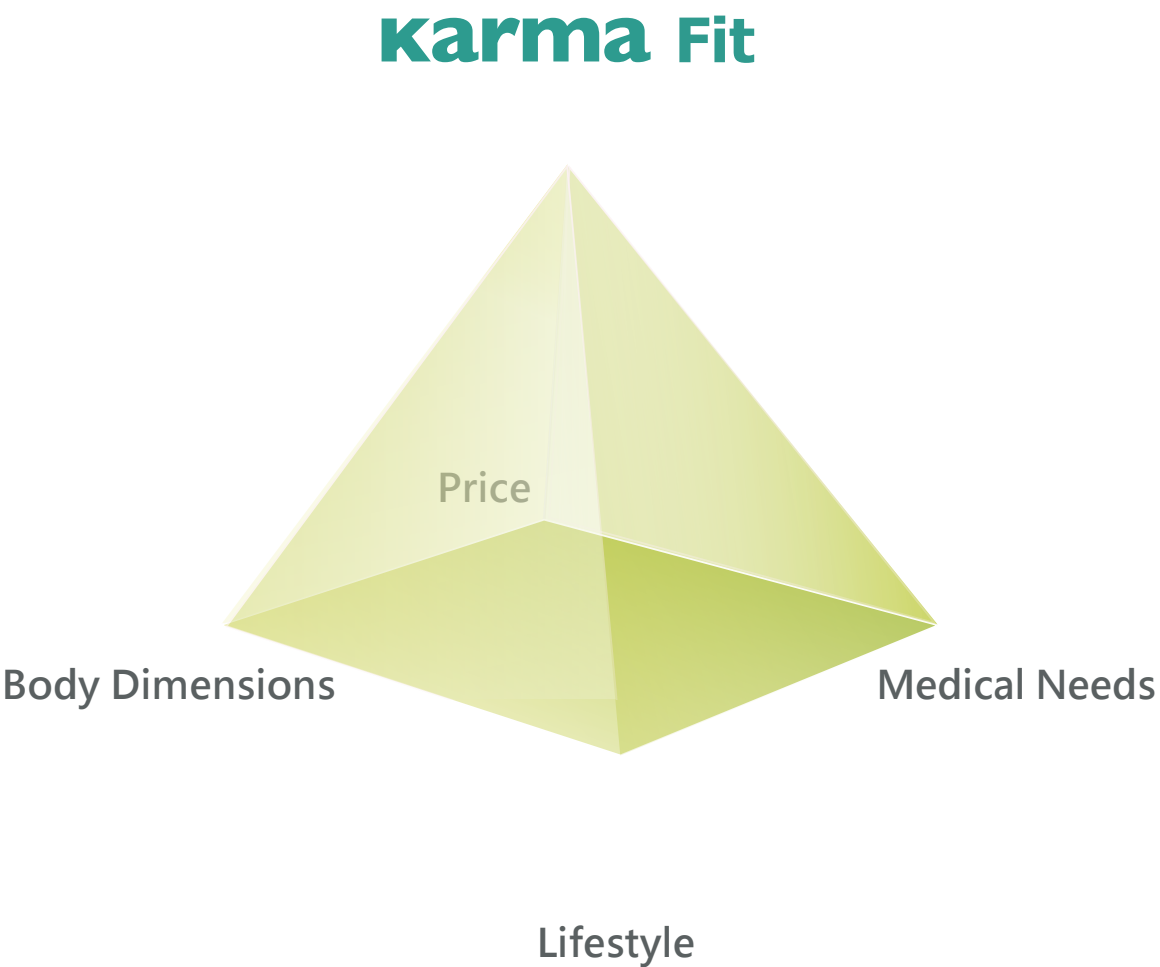

INDEX

S-Ergo Series

02
Series Intro

04
S-Ergo 125

05
S-Ergo 115

06
S-Ergo 105

07
S-Ergo 106

Ergo Lite Series

08
Series Intro

10
Ergo Lite

11
Ergo Lite 2

Ergo Semi-Active

12
Ergo Live

Flexx Series

16
Series Intro

20
Flexx / Flexx TALL

21
Flexx HD / Flexx S

Tilt-in-Space / Recline Series

22
Series Intro

23
VIP 515

24
MVP502

25
KM-5000

Heavy Duty

26
KM-BT10

Airport

27
KM-AA20

Pediatric Series

28
KM-7520

29
KM-CP33

Active Series

30
KM-AT60

Standard Series

32
Series Intro

33
KM-1500

34
KM-2500/KM-2500L

35
KM-TV20.2

Entry Level Series

36
Series Intro

37
SOMA 100/SOMA 105

38
SM-100.3/SM-150.3

39
ECON 800/ECON 805

Accessories


40

Specifiction

46

PRODUCT MAP

karma Manual Wheelchair

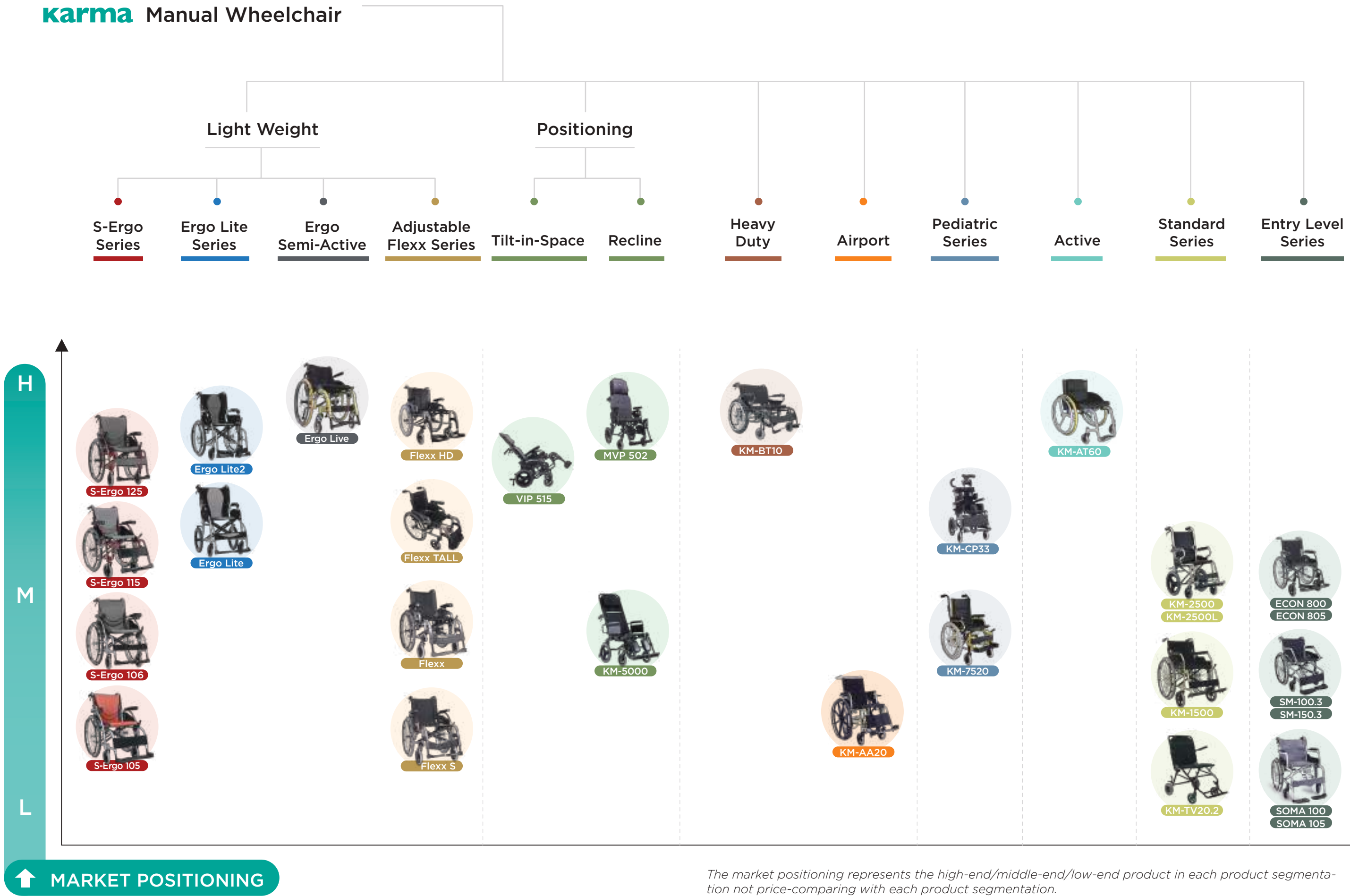

The market positioning represents the high-end/middle-end/low-end product in each product segmentation not price-comparing with each product segmentation.

S-Ergo Series

The patented S-Shaped Seating System comes with every wheelchair model in this series. We exclusively bring you a great product feature that can have a lasting effect on your everyday life.



Features for S-Ergo series

Shaped to fit you

The revolutionary S-Ergo system is the first s-shaped ergonomic seating in the world, designed according to the body curve, which is capable of simultaneously relieving pressure, stabilizing seating position and reducing downward sliding.






HOW DOES IT WORK?




TRADITIONAL SEAT

High risk of decubitus ulcers could be occurred because pressure points are concentrated on hipbone area.

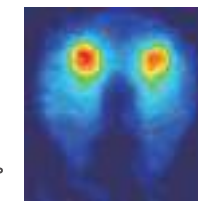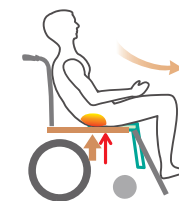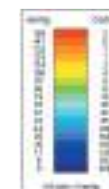

S-Ergo SEAT

The S-Ergo system provides efficient pressure relief by spreading weight over a greater area. This also increases stabilization and reduces sliding.

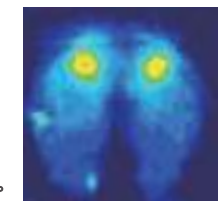

Ergo Armpad 2

The smooth edge of armpad provides comfortable resting and offer better hand movements during propelling.



Ergo Handrim

Patented ergonomic hand rim, easy to grip, easy to propel, and prevent from discomfort caused by the friction.



AEGIS

Anti-bacterial upholstery prevents from deterioration caused by bacteria, fungus and other microorganisms.



Compact

Light weight & compact, easy to be folded, convenient for storage and transportation.



Comfortable Seating and Easy Transfer

S-Ergo 125 is designed with detachable footrest and flip-back armrest, allowing users transfer easily in limited space. The innovative S-Ergo seating system provides comfortable seating environment.

**C.E.R.A.H**

Flip- back armrest offers greater convenience when the user transfer himself from wheelchair to bed, carseat or sofa.



Swing in/out and detachable footrest adds convenience when the user transfers in tight areas.



Quick release rear wheels reduce the weight and size for easy carrying and storage.

OPTIONS & ACCESSORIES

- Elevating Legrest
- Amputee Footrest
- Anti-tippers
- Pelvic Belt
- Harness
- Rigidifying Push Bar
- Foldable Headrest
- Lap Tray
- Travel Bag
- Backrest Bag
- Net Bag
- Cane Holder
- I.V. Pole Holder
- Canopy
- Cup Holder

FRAME COLOR

Rose Red Pearl Silver

SPECIFICATION

Seat Width (cm)	40 / 45 / 50
Seat Depth (cm)	43
Caster Size (inch)	7
Rear Wheel Size (inch)	14F / 24F / 24Q
Seat Height (cm)	50
Backrest Height (cm)	41
Armrest Height (cm)	21
Calf Length (cm)	31-40
Max. User Weight (kg)	115
Weight w/o Accessories (kg)	8.9
Total Weight (kg)	14.3
Overall L x W x H (cm)	99*58*91
Overall Collapsed L x W x H (cm)	80*31*74

±1cm ±0.5kg

-Karma reserves the right to modify information herein without further notice.

-For more details, please refer to the Specification Pages No. 46-53.

-F:fixed rear wheel axle / Q:quick release rear wheel axle



Stylish, Convenience, and Comfortable

S-Ergo 115 comes with 24" self propelling wheels and detachable footrest. Swing in/out and detachable footrest adds convenience when the user transfers in tight areas.

**C.E.R.A.H**

Swing in/out and detachable footrest adds convenience when the user transfers in tight areas.



Dual-purpose backrest/seat upholstery.



Foldable and extremely compact.

OPTIONS & ACCESSORIES

- Elevating Legrest
- Amputee Footrest
- Anti-tippers
- Pelvic Belt
- Harness
- Rigidifying Push Bar
- Foldable Headrest
- Lap Tray
- Travel Bag
- Backrest Bag
- Net Bag
- Cane Holder
- I.V. Pole Holder
- Canopy
- Cup Holder

FRAME COLOR

Rose Red Pearl Silver

SPECIFICATION

Seat Width (cm)	40 / 45 / 50
Seat Depth (cm)	43
Caster Size (inch)	7
Rear Wheel Size (inch)	14F / 24F / 24Q
Seat Height (cm)	50
Backrest Height (cm)	41
Armrest Height (cm)	21
Calf Length (cm)	31-40
Max. User Weight (kg)	115
Weight w/o Accessories (kg)	11.7
Total Weight (kg)	13.5
Overall L x W x H (cm)	99*59*91
Overall Collapsed L x W x H (cm)	80*31*74

±1cm ±0.5kg

-Karma reserves the right to modify information herein without further notice.

-For more details, please refer to the Specification Pages No. 46-53.

-F:fixed rear wheel axle / Q:quick release rear wheel axle

S-Ergo 105

C.E.R.A.H

Enjoy Unique Comfortable Seating

S-Ergo 105 is the standard version of S-Ergo series. It comes with effortless ergonomic brake and detachable seat covers. The streamlined oval tubes not only make the chair look great but also enhance the structure.



The outrigger armrest is designed to provide a comfortable arm resting support as well as maximize seating area.



14" rear wheel option

OPTIONS & ACCESSORIES

- Anti-tippers
- Pelvic Belt
- Harness
- Rigidifying Push Bar
- Foldable Headrest
- Lap Tray
- Backrest Bag
- Net Bag
- Cane Holder
- I.V. Pole Holder
- Canopy
- Cup Holder

FRAME COLOR

Rose Red Pearl Silver

SPECIFICATION

Seat Width (cm)	40 / 45
Seat Depth (cm)	43
Caster Size (inch)	7
Rear Wheel Size (inch)	14F / 24F / 24Q
Seat Height (cm)	50
Backrest Height (cm)	41
Armrest Height (cm)	21
Calf Length (cm)	31-40
Max. User Weight (kg)	115
Weight w/o Accessories (kg)	13.2
Total Weight (kg)	13.2
Overall L x W x H (cm)	102*58*91
Overall Collapsed L x W x H (cm)	98*31*74


±1cm ±0.5kg

-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 46-53.
-F:fixed rear wheel axle / Q:quick release rear wheel axle

S-Ergo 106

User-friendly, Comfortable and Easy to Steer

S-Ergo 106 adds extra support to the back. The easy adjustable backrest angle allows the user to cope with daily activities comfortably.

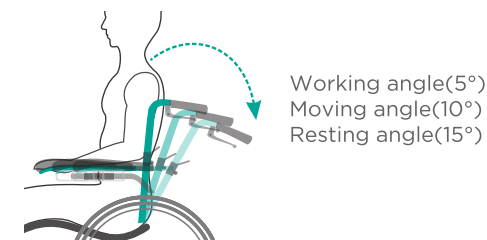

Tool-free angle adjustments for switching backrest between degree 5, 10, 15. The backrest-armrest interlinkage design allows the armrest self-adjusted while backrest angle changed.



Rigidify push bar strengthens back support, makes it easy to push and steer.



Anti-tippers ensure the safety of users.



OPTIONS & ACCESSORIES

- Pelvic Belt
- Harness
- Foldable Headrest
- Lap Tray
- Net Bag
- Cane Holder
- I.V. Pole Holder
- Canopy
- Cup Holder

FRAME COLOR

Pearl Silver

SPECIFICATION

Seat Width (cm)	45
Seat Depth (cm)	42
Caster Size (inch)	7
Rear Wheel Size (inch)	24F/24Q
Seat Height (cm)	50
Backrest Height (cm)	43
Armrest Height (cm)	23
Calf Length (cm)	31-40
Backrest Angle (X°)	5 / 10 / 15
Max. User Weight (kg)	115
Weight w/o Accessories (kg)	11
Total Weight (kg)	14.5
Overall L x W x H (cm)	104*63*90
Overall Collapsed L x W x H (cm)	98*31*74

±1cm ±0.5kg

-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 46-53.
-F:fixed rear wheel axle / Q:quick release rear wheel axle



Ergo Lite Series

For your on-the-go life!

Choose Ergo Lite ultra lightweight wheelchair series if you are looking for a portable option, folding frames collapse for everyday transport. It carries with innovative P&P Brake and patented S-Ergo Seating System, your best companion for travel.



Features for Ergo Lite series



S-Ergo seating system

Worldwide patented S-Ergo seating system relieves pressure of the bottom, reduces sliding and promotes good posture to prevent bedsores.



14" Rear Wheel Only

P&P brake system

Allow the caregiver to keep wheelchair standstill without bending over, increasing safety and convenience while going downhill.

Desk Armrest

Allow the user to get closer to the desk-top.



AEGIS

Anti-bacterial upholstery prevents from deterioration caused by bacteria, fungus and other microorganisms.



Reliable ultra lightweight concept

10kg is a fine line. It is said that the appropriate load for a woman to hand carry comfortably and with a lower risk of back injury is 10 kg. We strive to make the caretakers life easy. The Ergo Lite weighs only 8.5 kg complete and has successfully crash tested with capacity up to 100kg. The Ergo Lite 2 has even made the transportable piece to be 5.6kg.



Ultra Lightweight for Easy Travel

Ergo Lite is an ultra lightweight wheelchair that weight only 8.5kg. This compact design is ideal for female attendants. It is also equipped with Karma's unique S-Ergo seating system. We take care of both caregiver and user's need!



Only 8.5 kg, ultra-light weight and extremely compact.

Weight capacity 100 kg certified by CE regulations and UK Crash Test.



Convenient and suitable for caregivers to take care of elderly users. Easy to carry out and foldable for storage.

OPTIONS & ACCESSORIES

- Anti-tippers
- Pelvic Belt
- Harness
- Foldable Headrest
- Travel Bag
- Backrest Bag
- Cup Holder

SPECIFICATION

Seat Width (cm)	40 / 45
Seat Depth (cm)	42
Caster Size (inch)	6
Rear Wheel Size (inch)	14F
Seat Height (cm)	48
Backrest Height (cm)	42
Armrest Height (cm)	21
Calf Length (cm)	35-40
Max. User Weight (kg)	100
Weight w/o Accessories (kg)	8.5
Total Weight (kg)	8.5
Overall L x W x H (cm)	94*57*91
Overall Collapsed L x W x H (cm)	86*27*75

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 46-53.
-F:fixed rear wheel axle / Q:quick release rear wheel axle

FRAME COLOR

Pearl Silver

Best Companion for Lite Travel

Ergo Lite 2 has basic functions of the Ergo Lite but is improved with 20" rear wheels option and detachable footrest.



Self Propel



Transit



The weight without rear wheels and footrests is only 5.6kg, which is nearly equal to the weight of three hardcovers.



With swing-away and detachable footrests, helps users easily transfer out of the chair.



Optional 20" rear wheel with Ergo Handrim for self-propelling.

OPTIONS & ACCESSORIES

- Anti-tippers
- Pelvic Belt
- Harness
- Foldable Headrest
- Travel Bag
- Backrest Bag
- Net Bag
- Cup Holder

SPECIFICATION

Seat Width (cm)	40 / 45
Seat Depth (cm)	42
Caster Size (inch)	6
Rear Wheel Size (inch)	14F / 20F / 20Q
Seat Height (cm)	48
Backrest Height (cm)	42
Armrest Height (cm)	21
Calf Length (cm)	35-40
Max. User Weight (kg)	100
Weight w/o Accessories (kg)	7.7
Total Weight (kg)	9.2
Overall L x W x H (cm)	93*57*91
Overall Collapsed L x W x H (cm)	68*27*75

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 46-53.
-F:fixed rear wheel axle / Q:quick release rear wheel axle

FRAME COLOR

Pearl Silver



Ergo Live features crucial adjustments for personal comfort and fitness

Comfortable Seating



1 S-Ergo seating system provides pressure relief, stabilization and reduced sliding.



2 Unique design of backrest upholstery gives you extra support, also it provides AEGIS Microbe Shield.



3 Tension adjustable backrest for better back support.



4 Footrest angles and adjustable footplate allow the legs to be positioned comfortably.

Personal Configuration



1 Adjustable Rear Axle Position
Standard/ Prudential/ Active



2 Adjustable Backrest Angle
3-stage adjustment -5° /0° /5°



3 Adjustable Backrest Height
13"-16"/ 16"-19"



4 Adjustable Armrest Height



5 Adjustable Seat Height
17.5"-20"(with different
caster forks)



6 Multiple Seat Dimensions

Ergo Live

Ride with Style!

An enjoyable riding experience comes from a seamless integration of artistic design and delicate craftsmanship. Ergo Live's minimalist frames have fewer moving parts for maximum efficiency, especially designed for people who have great taste and wish to demonstrate personal style.



Options



5” Aluminum Rim Casters



24” Lively Rear Wheel (24-spoke)



24” Magnesium Rear Wheel



Push Handle



Foldable Push Handle



Skirt-mud Guard



Flip-back Armrest



2D Flip-up One Piece Footplate

OPTIONS & ACCESSORIES

- 5” Aluminum Rim Casters
- 24” Lively Rear Wheel (24-spoke)
- 24” Magnesium Rear Wheel
- Push Handle
- Foldable Push Handle
- Skirt-mud Guard
- Flip-back Armrest
- 2D Flip-up One Piece Footplate
- Anti-tippers
- Pelvic Belt
- Cane Holder

FRAME COLOR

- Lemon Yellow
- Fashion White
- Diamond Black

SPECIFICATION

Seat Width (cm)	36 / 41 / 46
Seat Depth (cm)	38 / 41 / 46
Caster Size (inch)	5 / 6
Rear Wheel Size (inch)	24Q
Seat Height (cm)	45-51
Backrest Height (cm)	33-40.5/40.5-48
Armrest Height (cm)	23
Calf Length (cm)	36-44
Backrest Angle(X°)	(-5) / 0 / 5
Max. User Weight (kg)	125
Weight w/o Accessories (kg)	8.9
Total Weight (kg)	11.9
Overall L x W x H (cm)	Min. 89*56*83 / Max. 89*56*97
Overall Collapsed L x W x H (cm)	Min. 89*32*83 / Max. 89*32*97

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 46-53.
-F:fixed rear wheel axle / Q:quick release rear wheel axle





Flexx Series

Robust, Stylish and Adjustable.



Choose Flexx series if you are an active user looking for an ideal custom-made wheelchair. Flexx is multi-adjustable to fulfill your need, as an extension of you, it is highly responsive, and offer precise control, superb balance, and exceptional performance.



Features for Flexx Series

Configurable and Adjustable

Ideal for institutional center usage, which can be adjusted for various users.

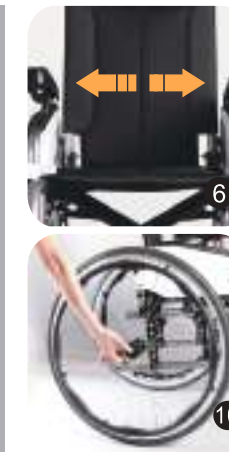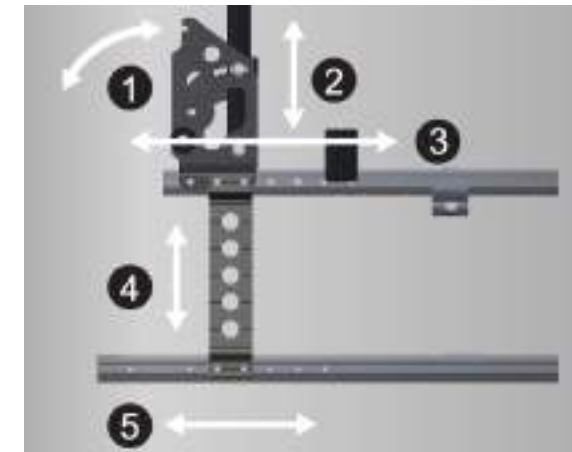

- 1 Backrest angles
- 2 Backrest heights
- 3 Seat depths
- 4 Seat to ground heights
- 5 Adjustable axle plate
- 6 Seat widths
- 7 Front caster
- 8 Rear wheel
- 9 25 Rear wheel positions
- 10 Quick release axle

Universal Design

The bright color of key lockers is eye-catching and understandable, users can easily operate the functions without second thoughts.



Transfer with ease

Flip-back armrest and swing-away footrest improve convenience and easier access to the chair.



Patented Ergonomic Hand Rim

Patented ergonomic hand rim, easy to grip, easy to propel, and prevent from discomfort caused by the friction.



AEGIS

Anti-bacterial upholstery prevents from deterioration caused by bacteria, fungus and other microorganisms.

Flexx provides versatile options to satisfy individual needs

Armrest

Ergo Armpad 2 (Standard)

Boat Armpad

Hemi Armrest

Legrests & Footplates

Elevating Legrest
Offer optimal compensation with adjustable lengths from 38-48 cm.

Amputee Legrest
Height, angle and depth adjustable for perfect support.

Aluminium Footplate

Legrests & Footplates

90° Legrest & Multi-directional Adjustable Footplate
The functions of 90° Legrests and depth / angle adjustable footplate can match different users' leg shape comfortably.

Angle Adjustable Footplate
70°~110° angle adjustable

1" & 2" Extended Tube
Seat depth can be added 1" or 2" for ideal sitting comfort without replacement of seat cushion.

Backrest & Seat

Foldable Backrest

Adjustable Backrest Height
Backrest height can be adjusted 38-53cm

Tension Adjustable Backrest Upholstery
Prevent a hammock effect. Can be shaped to enhance comfort and fit individual users needs.

Fork

Short/Long Fork (Standard)

Stylish Short/Long Fork

2" Wide Fork

Caster

5"

6"

7"

8"

200x50

Headrest

Foldable Headrest

Rear Wheel

14"

20"

22"

24"

Quick Release Wheels

Safety Options

Tie-Down Bracket

Swing-away&Height Adjustable Anti Tipper
Easily folded away using the foot without bending over.

Brake

Attendant Brake



Fully adjustable for the best comfort and propelling performance.



Extended seat depth to assist taller users in need of mobility.



Extra durable and upgraded capacity for plus-size individuals in need.



The entry level of Flexx series, but equipped with all key adjustments and features.



OPTIONS & ACCESSORIES	FLEXX	FLEXX TALL	FLEXX HD	FLEXX S
• Foldable Backrest	○	○	○	○
• Stylish Fork	○	○	●	-
• Tension Adjustable Backrest	○	○	○	○
• 1 1/2" Extended Seat Tube	○	-	-	○
• Boat Armpad	○	○	○	○
• Hemi Armpad	○	○	○	○
• Axle Traveling Elevating Legrest	○	○	○	○
• Amputee Footrest	○	○	○	○
• Vertical Footrest	○	-	-	○
• Swing-away Anti-tippers	○	○	○	○
• Attendant Brake	○	○	○	○
• Tie-down Bracket	○	○	○	○
• Pelvic Belt	○	○	○	○
• Harness	○	○	-	○
• Rigidifying Push Bar	○	○	●	○
• Foldable Headrest	○	○	-	○
• Gel Cushion	○	-	-	○
• Lap Tray	○	○	-	○
• Backrest bag	○	○	-	○
• Net bag	○	○	-	○
• Cane Holder	○	○	○	○
• I.V. Pole Holder	○	○	○	○
• Oxygen Tank Holder	○	○	○	○
• Cup Holder	-	-	-	○

● Standard ○ Option - Not available

SPECIFICATION

	FLEXX	FLEXX TALL	FLEXX HD	FLEXX S
Seat Width (cm)	38 / 40.5 / 43 / 45.5 / 48 / 50.5	40.5/43/45.5/48	51 / 56 / 61	38 / 40.5 / 43 / 45.5 / 48 / 50.5
Seat Depth (cm)	41 / 43.5 / 46	48.5/51	41 / 46 / 51	41 / 43.5 / 46
Caster Size (inch)	5 / 6 / 7 / 8	6 / 7 / 8	6 / 7 / 8	5 / 6 / 7 / 8
Rear Wheel Size (inch)	14/20/22/24	14/22/24	14/22/24	14/20/22/24
Seat Height (cm)	39/40.5/43/45.5/48/50.5/53	45.5/48/50.5/53	45.5/48/50.5/53	39/40.5/43/45.5/48/50.5/53
Backrest Height (cm)	43/51/38-53	43/51/38-46	43/51/38-46	43/51/38-53
Armrest Height (cm)	21-29	21-29	21-29	22-24
Calf Length (cm)	38-48	38-48	38-48	38-48
Backrest Angle(X°)	(-4)-16	(-4)-16	(-4)-16	(-4)-16
Seat Angle (X°)	3-7	3-7	3-7	3-7
Max. User Weight (kg)	130	130	170	130
Weight w/o Accessories (kg)	Min. 9.3	Min. 10.6	Min. 10.6	Min. 9.3
Total Weight (kg)	Min. 14.2 (14F) / Min. 14.7 (24Q)	Min. 15.2 (14F) / Min. 15.7 (24Q)	Min. 16.5 (14F) / Min. 17.0 (24Q)	Min. 14.2 (14F) / Min. 14.7 (24Q)
Overall L x W x H (cm)	Min. 102*56*76 / Max. 107*69*106	Min. 115*59*84 / Max. 118*67*104	Min. 108*71*84 / Max. 118*81*104	Min. 102*56*76 / Max. 107*69*106
Overall Collapsed Width (cm)	34	34	35	34

±1cm ±0.5kg

-Karma reserves the right to modify information herein without further notice.

-For more details, please refer to the Specification Pages No. 46-53.

-F:fixed rear wheel axle / Q:quick release rear wheel axle


FRAME COLOR Silver Black

Tilt-in-Space / Recline Series

Discover your relaxed manner!

Tilt-in-space/Recline series wheelchairs are designed for a variety of positioning needs. The tilted position helps to facilitate feeding and respiratory function, reduce pressure beneath the pelvis and improve visual alignment by holding the head upright.

Release your body pressure

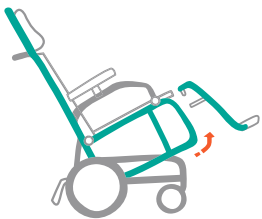

Tilt-in-Space | VIP 515

For users without self-supporting capacity


Tilt-in-Space function effectively release bottom pressure, easier and safer to improve user's positioning.



▲ Pressure is centralized on normal wheelchairs




▲ Pressure is distributed when using tilt-in-space

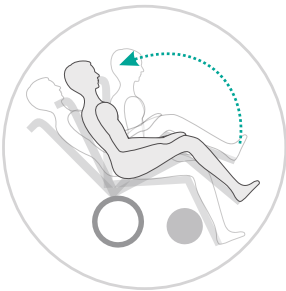

▲ Center-of-gravity weight shifting tilt-in-space

Intelligent center-of-gravity weight shifting system






This unique system allows center of gravity to shift between front and rear wheels while tilted up, providing the steadiest structure possible, along with greater safety as well.

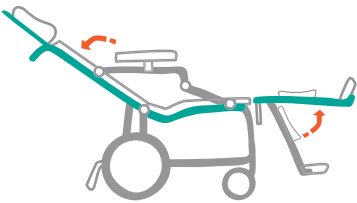

Recline | MVP 502 /KM-5000

For users with limited supporting capacity

With Reclining function, caregiver is able to proceed backslap, massage and diaper change with ease.

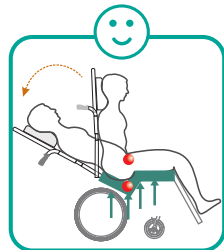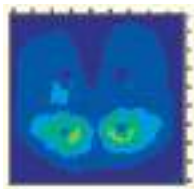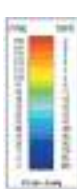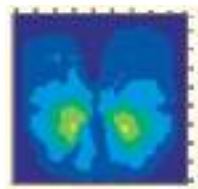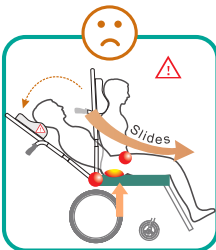

Patented anti-sliding V-seating system (MVP502)

The V-shaped seat smoothly shifts bottom 20° downward, minizing shear force when reclining. It not only prevents sliding, but also reduces the risk of decubitus ulcers.

VIP 515



Redefining the Concept of Comfort

VIP 515 is a revolutionary lightweight folding tilt-in-space wheelchair, with a tilt range from 0 to 35 degrees. Carried exclusive gas struts for an intelligent center-of-gravity weight shifting system to maintain stability while the chair is tilted up.



Compact, easy to fold, store and transport.



Flip-back & height adjustable armrests, swing-away & detachable footrests provides easier access to the chair.



20" rear wheel option



Ergonomic headrest relieves tension on shoulder and neck.

OPTIONS & ACCESSORIES

- 1"/2" Extended Seat Tube
- Boat Armpad
- Axle traveling Elevating Legrest
- Amputee Footrest
- Vertical Footrest
- Tie-down Bracket
- Harness
- Rigidifying Push Bar
- Gel Cushion
- Lap Tray
- Backrest Bag
- Net bag
- Cane Holder
- I.V. Pole Holder
- Oxygen Tank Holder
- Cup Holder

FRAME COLOR

Diamond Black

SPECIFICATION

Seat Width (cm)	40 / 45
Seat Depth (cm)	42
Caster Size (inch)	7
Rear Wheel Size (inch)	14F / 20F / 20Q
Seat Height (cm)	46
Backrest Height (cm)	49
Armrest Height (cm)	21-29
Calf Length (cm)	41-46
Seat Angle (X°)	0-35
Max. User Weight (kg)	115
Weight w/o Accessories (kg)	15.7
Total Weight (kg)	18.9
Overall L x W x H (cm)	108*62*99(134)
Overall Collapsed L x W x H (cm)	73*34*76

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 46-53.
-F:fixed rear wheel axle / Q:quick release rear wheel axle

MVP 502




Designed for Comfort and Posture Stabilization

MVP 502 is with innovated anti-sliding V-shaped seating, it effectively releases pressure for bottom, greatly reduces the unpleasant feeling from back after position change.



Reclining angle from 0 to 67 degrees.

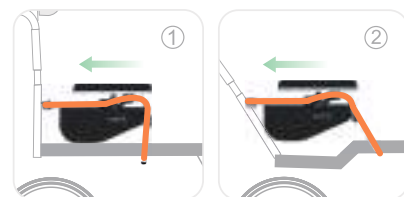

The backrest-armrest coupling structural design.



Axle-traveling elevating footrests make the legs to be stretched smoothly without compression.



Easy to fold, store and transport.

OPTIONS&ACCESSORIES

- 1"/2" Extended Seat Tube
- Amputee Footrest
- Vertical Footrest
- Boat Armpad (For 14" rear wheel only)
- Harness
- Lap Tray
- Net bag
- Cane Holder
- I.V. Pole Holder
- Oxygen Tank Holder

FRAME COLOR

Diamond Black

SPECIFICATION

Seat Width (cm)	37 / 40 / 42 / 45 / 50
Seat Depth (cm)	42 / 46
Caster Size (inch)	7
Rear Wheel Size (inch)	14F / 22F / 22Q
Seat Height (cm)	51
Backrest Height (cm)	51
Armrest Height (cm)	23-29
Calf Length (cm)	42-49
Backrest Angle(X°)	0-67
Legrest Angle (X°)	0-70
Max. User Weight (kg)	115
Weight w/o Accessories (kg)	15
Total Weight (kg)	20
Overall L x W x H (cm)	118*57*96(134)
Overall Collapsed L x W x H (cm)	85*30*76

±1cm ±0.5kg

-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 46-53.
-F:fixed rear wheel axle / Q:quick release rear wheel axle

KM-5000



Standard Version of Recline Series

KM-5000 is a standard version of recline series, especially designed for long-term user and stroke patients, it's full detachable armrests and elevating leg rest offer great convenience to transfer, ergonomic headrest and adjustable backrest allows users to sit and lie down comfortably.



Backrest can be reclined smoothly from 5 to 75 degrees.



Full detachable armrests offer great convenience to transfer from the wheelchair.



With elevating & swing-away footrests.



20" rear wheel option

OPTIONS&ACCESSORIES

- Amputee Footrest
- Vertical Footrest
- Attendant Brake
- Pelvic Belt
- Harness
- Gel Cushion
- Lap Tray
- Backrest Bag
- Net bag
- Cane Holder
- I.V. Pole Holder
- Oxygen Tank Holder

FRAME COLOR

Black

SPECIFICATION

Seat Width (cm)	41 / 46 / 51
Seat Depth (cm)	43 / 45
Caster Size (inch)	8
Rear Wheel Size (inch)	14F / 24F / 24Q
Seat Height (cm)	51
Backrest Height (cm)	45
Armrest Height (cm)	22
Calf Length (cm)	41-48
Backrest Angle (X°)	5-75
Legrest Angle (X°)	0-75
Max. User Weight (kg)	115
Weight w/o Accessories (kg)	15
Total Weight (kg)	20.6
Overall L x W x H (cm)	127*58*93(129)
Overall Collapsed L x W x H (cm)	91*34*93

±1cm ±0.5kg


-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 46-53.
-F:fixed rear wheel axle / Q:quick release rear wheel axle

KM-BT10

Weight Capacity up to 250 kg with Highly Adjustable and Foldable Design

KM-BT10 is with sturdy but foldable frame, exclusively supporting heavy or plus-size individuals in need, capacity up to 250kg, max. seat width 76cm.

Prescription based seat and back system.

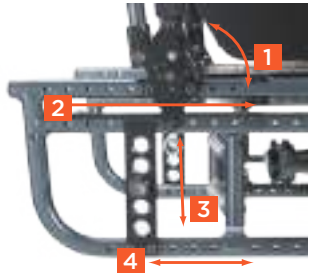


- 1 Back Angle Adjustment
- 2 Seat depth
- 3 Seat to Floor Height
(Available with different wheel and caster package)
- 4 Adjustable axle plate



Solid seat pan offers great support and prevent Hammock Effect.



Unique cross brace provides robust frame structure and makes it also foldable.

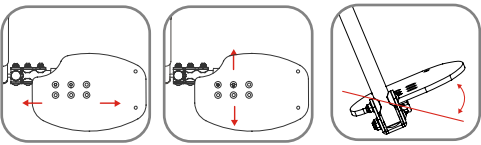

Multi-directional adjustable footplates.

OPTIONS & ACCESSORIES

- Height Adjustable Armrest
- Elevating Legrest
- Adjustable Heel Loop
- 7.5 cm Cushion
- Quick release rear wheel
- 5mm-space Handrim
- I.V. Pole Holder
- Oxygen Tank Holder

FRAME COLOR

Diamond Black

SPECIFICATION

Seat Width (cm)	56 / 61 / 66 / 71 / 76
Seat Depth (cm)	46 - 51 / 56 - 61
Caster Size (inch)	7 / 8
Rear Wheel Size (inch)	20F/20Q/24F/24Q
Seat Height (cm)	39-50
Backrest Height (cm)	41-46/46-51
Armrest Height (cm)	25
Calf Length (cm)	32-37
Backrest Angle (X°)	(-5)-15
Max. User Weight (kg)	250
Weight w/o Accessories (kg)	22.9
Total Weight (kg)	31
Overall L x W x H (cm)	Min. 118*78*80 Max. 123*88*101
Overall Collapsed L x W x H (cm)	Min. 90*44*80 Max. 95*44*101

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 46-53.
-F:fixed rear wheel axle / Q:quick release rear wheel axle

KM-AA20

Easy to Travel in Airport or Narrow Space

KM-AA20 is designed for airplane transport, ideal for use in airport environment, such as terminal lobby, corridor. After rear wheel is released, it can travel freely inside the aircraft, even in narrow aisle.



Sturdy non-foldable frame. With quick release rear wheel, users can easily shift between big and small rear wheels and without any tools required.

Only 54cm wide once the rear wheels are detached.



Swing-away armrests provide easily lateral transfer.

OPTIONS & ACCESSORIES

- Harness

FRAME COLOR

Champagne

SPECIFICATION

Seat Width (cm)	39
Seat Depth (cm)	42
Caster Size (inch)	6
Rear Wheel Size (inch)	24Q
Seat Height (cm)	49
Backrest Height (cm)	42
Armrest Height (cm)	19
Calf Length (cm)	32
Max. User Weight (kg)	100
Weight w/o Accessories (kg)	11.5
Total Weight (kg)	16.2
Overall L x W x H (cm)	99*58*91
Overall Collapsed L x W x H (cm)	99*58*91

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 46-53.
-F:fixed rear wheel axle / Q:quick release rear wheel axle

KM-7520

Multi-functional to Fit Children's Needs

KM-7520 is featured with adjustable design, supporting the user's body, life style, environment and skill changes.



Height adjustable push handles make it push and steer easily.



With full detachable armrest and swing-away footrest.



Height adjustable armrests provide optimal comfort for child at any moment.



Optional 20" quick release rear wheels meet the demand of both self-propelling and minimized storage space.



Recline

Optional reclining version to help bottom pressure relief.

OPTIONS & ACCESSORIES

- Elevating Legrest
- Pelvic Belt
- Lap Tray
- Cane Holder
- I.V. Pole Holder
- Oxygen Tank Holder

SPECIFICATION

	KM-7520	KM-7520 Recline
Seat Width (cm)	28 / 30.5 / 34	28 / 30.5 / 34
Seat Depth (cm)	28 / 30.5 / 34	28 / 30.5 / 34
Caster Size (inch)	6	6
Rear Wheel Size (inch)	12F / 20Q	20Q
Seat Height (cm)	46	46
Backrest Height (cm)	41	41
Armrest Height (cm)	17-23	17-23
Calf Length (cm)	29	29
Backrest Angle (X°)	-	0-60
Legrest Angle (X°)	-	0-75
Max. User Weight (kg)	100	100
Weight w/o Accessories (kg)	12	13.5
Total Weight (kg)	14.5	17
Overall L x W x H (cm)	81*51*91	85*51*91
Overall Collapsed L x W x H (cm)	62*33*84	62*33*84

FRAME COLOR

Yellow

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 46-53.
-F:fixed rear wheel axle / Q:quick release rear wheel axle

KM-CP33

Perfect Pressure-free Sitting Experience for Children

The combination of Tilt-in-space and Recline system on KM-CP33 provides a perfect pressure-free sitting environment for children. Separator, trunk supporter, dining tray and many other accessories and options are available to personalize the most suitable KM-CP33 for each user.



Lightweight frame is easy to store and transport.



Customized design fits users' body condition and medical needs.



Tilt-in-Space system distributes bottom pressure, stabilizes the sitting posture.



Recline system allows kids to stretch abdominal muscle and reduces the stress and discomfort.

OPTIONS & ACCESSORIES

- Foam Headrest
- Claw Shape Headrest
- Curved Shape Headrest
- H-shaped Chest Belt
- Pelvic Belt
- Trunk Support
- Thigh Support
- Thigh Separator (Pommel)
- Arm Loops/Thigh Loops/Calf Loops/Heel Loops/Thenar Loops
- Lap Tray
- I.V. Pole Holder
- Oxygen Tank Holder

SPECIFICATION

Seat Width (cm)	20.5/ 23/ 25.5/ 28/ 30.5/ 33/ 35.5/ 38
Seat Depth (cm)	23/ 25.5/ 28/ 30.5/ 33/ 35.5/ 38/ 40.5
Caster Size (inch)	6
Rear Wheel Size (inch)	14F / 20F
Seat Height (cm)	45
Backrest Height (cm)	25-56
Armrest Height (cm)	15-33
Calf Length (cm)	38-46
Backrest Angle (X°)	0-30
Seat Angle (X°)	0-40
Legrest Angle (X°)	0-75
Max. User Weight (kg)	80
Weight w/o Accessories (kg)	Min. 9
Total Weight (kg)	Min. 11
Overall L x W x H (cm)	Min. 83*43*65 Max. 106*61*96
Overall Collapsed Width (cm)	40

FRAME COLOR

Black

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 46-53.
-F:fixed rear wheel axle / Q:quick release rear wheel axle



Active series


Explore your limits!

Active series wheelchairs are designed with ultra lightweight materials, high performance technology, and efficient ergonomics to bring you a competitive edge.

KM-AT60

Alluring, Stylish and Easy to Carry

With the superior agility and maneuverability design, the KM-AT60 is especially featured with ultra-light weight and extra strength. This model conquers the difficulties of non-foldable frames of active chair and makes it easy to carry and transport.

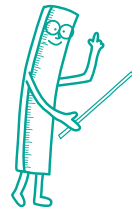

Customized design for seat width, seat height, backrest height, backrest angle, footrest height, camber and adjustable Center of Gravity. Various options for rear wheels and other accessories!



Rigid frame with foldable back & quick release rear wheel design facilitates the daily transportation and storage.



Tension adjustable backrest prevents Hammock Effect and gives users optimal back support.



Various options for rear wheels and other accessories.

OPTIONS & ACCESSORIES

- Sport Rear Wheel (Red wheel rim)
- Magnesium Rear Wheel (Silver or Black rim)
- Adjustable Push Handle
- Clip-on Push Handle
- Scissor Brake
- Anti-tippers
- Leg-tube Protector
- Spoke Guard
- 5 cm Cushion

SPECIFICATION

Seat Width (cm)	38 / 40.5 / 43 / 45.5
Seat Depth (cm)	35.5 / 38 / 40.5 / 43 / 45.5 / 48 / 50.5
Caster Size (inch)	5
Rear Wheel Size (inch)	24Q
Seat Height (cm)	46-58
Backrest Height (cm)	31-48
Calf Length (cm)	33-48
Seat Angle (X°)	0 / 3 / 5
Camber angle (X°)	0 / 3 / 5 / 7
Max. User Weight (kg)	100
Weight w/o Accessories (kg)	Min. 8
Total Weight (kg)	Min. 10.5
Overall L x W x H (cm)	Min. 88*62*77 Max. 103*70*106

FRAME COLOR

- Rose Red
- Swan Yellow
- Diamond Black
- Diamond Blue

±1cm ±0.5kg
 -Karma reserves the right to modify information herein without further notice.
 -For more details, please refer to the Specification Pages No. 46-53.
 -F:fixed rear wheel axle / Q:quick release rear wheel axle



Standard Series

Start from versatile standard

Karma standard series are made from durable aluminum folding frames. They are high quality, reliable, and include options for the comfort of users and caregivers alike.

KM-1500



Standard Version and Easy Storage

KM-1500 is an ideal choice in its category for standard body type and operating environment . Appropriate for general users, healthcare centers and hospitals. With a lightweight aluminum frame and foldable components, the KM-1500 helps to function at your best and responds to your mobility needs.



The frame is made of lightweight aircraft-grade aluminum alloy. It's a Karma classical wheelchair with fixed armrest and fixed footrest.



Foldable backrest reduces storage space and easy to carry out.



Outward swivel footplates provides more space for movement and transfer.



Dual brake system brings more convenience for caregivers and enhances the safety.

OPTIONS & ACCESSORIES

- Anti-tippers
- Pelvic Belt
- Harness
- Foldable Headrest
- Gel Cushion
- Lap Tray
- Backrest Bag
- Net Bag
- Cane Holder
- I.V. Pole Holder
- Oxygen Tank Holder
- Cup Holder

SPECIFICATION

Seat Width (cm)	39 / 45
Seat Depth (cm)	42
Caster Size (inch)	8
Rear Wheel Size (inch)	24F
Seat Height (cm)	48
Backrest Height (cm)	38
Armrest Height (cm)	21
Calf Length (cm)	32-43
Max. User Weight (kg)	100
Weight w/o Accessories (kg)	13.5
Total Weight (kg)	13.5
Overall L x W x H (cm)	100*63*87
Overall Collapsed L x W x H (cm)	91*31*72

FRAME COLOR

- Champagne
- Black

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 46-53.
-F:fixed rear wheel axle / Q:quick release rear wheel axle

KM-2500/ KM-2500L

10.9 kg Super Lightweight for Easy Transport
KM-2500/KM-2500L is an ideal choice in its category for people with thin, small body type or in narrow operating environment. Lightweight aluminum frame and foldable design is space saving and easy for travelling.



Aluminum alloy frame, weight only 10.9kgs.

 Compact and foldable backrest reduces storage space and easy for traveling.


 Outward swivel footplates provide more space for movement and transfer.

 Optional 20" rear wheel with Ergo Handrim for self-propelling. (KM-2500L)

OPTIONS & ACCESSORIES

- Anti-tippers
- Attendant Brake
- Pelvic Belt
- Harness
- Foldable Headrest
- Gel Cushion
- Lap Tray
- Travel Bag
- Backrest Bag
- Net Bag
- Cane Holder
- I.V. Pole Holder
- Cup Holder

FRAME COLOR

-  Champagne
-  Black

SPECIFICATION

	KM-2500	KM-2500L
Seat Width (cm)	41 / 47	41 / 47
Seat Depth (cm)	40	40
Caster Size (inch)	6	6
Rear Wheel Size (inch)	14F	20F / 20Q
Seat Height (cm)	47	47
Backrest Height (cm)	39	39
Armrest Height (cm)	21	21
Calf Length (cm)	32-43	32-43
Max. User Weight (kg)	100	100
Weight w/o Accessories (kg)	10.9	11.7
Total Weight (kg)	10.9	11.7
Overall L x W x H (cm)	89*53*86	89*56*86
Overall Collapsed L x W x H (cm)	72*25*68	82*28*68

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 46-53.
-F:fixed rear wheel axle / Q:quick release rear wheel axle


KM-TV20.2

Ideal Choice for Travelling
KM-TV20.2 is an excellent choice for getting the maximum convenience and comfort whether traveling to and from the house, to and from the car, or doing short errands.



As compact as a stack of large pizza boxes after folded, easy to carry out and enjoy a pleasant journey.

3 simple steps to fold the chair up and without any tools required.

 Padded flip-up armrests.

 Secure brake improves safety.

OPTIONS & ACCESSORIES

- Harness

SPECIFICATION

Seat Width (cm)	41/46
Seat Depth (cm)	47
Caster Size (inch)	6
Rear Wheel Size (inch)	8
Seat Height (cm)	53
Backrest Height (cm)	47
Armrest Height (cm)	28
Max. User Weight (kg)	100
Weight w/o Accessories (kg)	8.9
Total Weight (kg)	8.9
Overall L x W x H (cm)	91*60*92
Overall Collapsed L x W x H (cm)	61*35*71

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 46-53.
-F:fixed rear wheel axle / Q:quick release rear wheel axle

FRAME COLOR

-  Black



Entry Level Series

Entry level series are designed to make caregiver-assisted travel safe, convenient, and comfortable. Whether traveling to and from the house, to and from the car, or doing short errands, these wheelchairs offer comfortable sitting surfaces, narrow frames, and fold easily for storage.

SOMA 100/ SOMA 105

Simplicity but Not Simple

SOMA 100/SOMA105 are crafted from durable frames and are perfect for occasional or short-term use. Whether propelling the wheelchair using your arms, feet, or with assistance from a caregiver, SOMA wheelchairs include adjustable components and folding frames for comfort and convenience.



Modern, stylish and unique elliptical design with streamline shape.



Ergonomic push handle provides a stable and easy operation for attendant braking.



Energy-saving tipping lever helps caregivers to overcome obstacles easier.



16" rear wheel option

OPTIONS & ACCESSORIES

- Foldable Backrest (SOMA 105)
- Attendant Brake
- Pelvic Belt
- Harness
- Gel Cushion
- Lap Tray
- Backrest Bag
- Net Bag
- Cane Holder
- I.V. Pole Holder
- Cup Holder

SPECIFICATION

	SOMA 100	SOMA 105
Seat Width (cm)	40/45	40/45
Seat Depth (cm)	42	42
Caster Size (inch)	7	7
Rear Wheel Size (inch)	16F/22F	16F/22F
Seat Height (cm)	50	50
Backrest Height (cm)	38	38
Armrest Height (cm)	22	22
Calf Length (cm)	41-45	41-45
Max. User Weight (kg)	100	100
Weight w/o Accessories (kg)	11.3	11.7
Total Weight (kg)	11.3	11.7
Overall L x W x H (cm)	97*68*87	97*68*87
Overall Collapsed L x W x H (cm)	97*28*87	97*28*73

±1cm ±0.5kg
 -Karma reserves the right to modify information herein without further notice.
 -For more details, please refer to the Specification Pages No. 46-53.
 -F:fixed rear wheel axle / Q:quick release rear wheel axle

UPHOLSTERY COLOR

Black Check Red Check

FRAME COLOR

Silver White

SM-100.3/
SM-150.3

Standard Version and Fit Your Budget

SM-100.3/SM-150.3 are built using lightweight aluminum and are great options when you need to use a wheelchair more frequently, want a intermediate level of adjustability. Both two are crafted from durable frames and are perfect for occassional or short-term use.



Standard wheelchair but with unique and stylish elliptical frame structure.



Detachable seat upholstery design is convenient for cleaning.

Integrated seat system prevents Hammock Effect.



Energy-saving tipping lever helps caregivers to overcome obstacles easier.



Anti-slip and ergonomic footplates provide a stable and comfortable base for feet.

OPTIONS & ACCESSORIES

- Foldable Backrest (SM-150.3)
- Attendant Brake
- Pelvic Belt
- Harness
- Gel Cushion
- Lap Tray
- Backrest Bag
- Net Bag
- Cane Holder
- I.V. Pole Holder
- Cup Holder

SPECIFICATION

	SM-100.3	SM-150.3
Seat Width (cm)	43	43
Seat Depth (cm)	42	42
Caster Size (inch)	7	7
Rear Wheel Size (inch)	16F / 22F	16F / 22F
Seat Height (cm)	51	51
Backrest Height (cm)	38	38
Armrest Height (cm)	23	23
Calf Length (cm)	41-44	41-44
Max. User Weight (kg)	100	100
Weight w/o Accessories (kg)	11.3	11.5
Total Weight (kg)	11.3	11.5
Overall L x W x H (cm)	99*64*88	99*64*88
Overall Collapsed L x W x H (cm)	99*29*88	99*29*73

FRAME COLOR

- SM-100.3 Black
- SM-150.3 Silver White

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 46-53.
-F:fixed rear wheel axle / Q:quick release rear wheel axle

ECON 800/
ECON 805

Economical Version and Basic Transfer Function

ECON 800/ECON 805 are designed to make caregiver-assisted travel safe, convenient, and comfortable. Lightweight is easy to store and tranport, they are highly responsive, lightweight and offer precise control, superb balance and exceptional performance.



Aluminum alloy frame strikes a perfect balance between weight and strength.

Padded armrests for better comfort.



Flip-back armrests help users easily transfer from wheelchair to other places.



Detachable and swing-away footrests offer easy access to the chair.



Optional quick release 24" rear wheels to reduce storage space.

OPTIONS & ACCESSORIES

- Foldable Backrest (ECON 805)
- Elevating Legrest
- Amputee Footrest
- Anti-tippers
- Attendant Brake
- Pelvic Belt
- Harness
- Foldable Headrest
- Gel Cushion
- Lap Tray
- Backrest Bag
- Net Bag
- Cane Holder
- I.V. Pole Holder
- Oxygen Tank Holder

SPECIFICATION

	ECON 800	ECON 805
Seat Width (cm)	41/46/51	41/46/51
Seat Depth (cm)	43	43
Caster Size (inch)	8	8
Rear Wheel Size (inch)	16F/24F/24Q	16F/24F/24Q
Seat Height (cm)	48	48
Backrest Height (cm)	43	43
Armrest Height (cm)	23	23
Calf Length (cm)	39-48	39-48
Max. User Weight (kg)	115	115
Weight w/o Accessories (kg)	12.1	12.4
Total Weight (kg)	14.6	14.9
Overall L x W x H (cm)	111*59*94	111*59*94
Overall Collapsed L x W x H (cm)	87*32*94	87*32*75

FRAME COLOR

- Black

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 46-53.
-F:fixed rear wheel axle / Q:quick release rear wheel axle

Accessories

*The product pictures are for reference only, the real goods varied according to customer requirements.

Foldable Backrest



Standard For Models

S-Ergo 125 MVP 502
S-Ergo 115 KM-CP33
S-Ergo 105 KM-1500
S-Ergo 106 KM-2500/KM-2500L
Ergo Lite SOMA 105
Ergo Lite 2 SM-150.3
VIP 515 ECON 805

Optional For Models

Flexx S
Flexx
Flexx Tall
Flexx HD

Stylish Fork



Standard For Models

Ergo Live
Flexx HD

Optional For Models

Flexx
Flexx Tall

Tension Adj. Backrest Upholstery



Standard For Models

Ergo Live
MVP 502
KM-BT10
KM-AT60

Optional For Models

Flexx S
Flexx
Flexx Tall
Flexx HD

1"/2" Extended Seat Tube



Standard For Models

None

Optional For Models

Flexx S
Flexx
VIP 515 (with Footrest Hanger Extension Bracket)
MVP 502 (with Footrest Hanger Extension Bracket)

Footrest Hanger Extension Bracket



Standard For Models

None

Optional For Models

VIP 515
MVP 502

Flip-back and Height-adj. Armrest



Standard For Models

Flexx S(by tool)
Flexx
Flexx Tall
Flexx HD
VIP 515
MVP 502
KM-7520(Full Detachable)

Optional For Models

KM-BT10

Boat Armpad



Standard For Models

None

Optional For Models

Flexx S
Flexx
Flexx Tall
Flexx HD
VIP 515
MVP 502

Hemi Armpad



Standard For Models

None

Optional For Models

Flexx S
Flexx
Flexx Tall
Flexx HD

Elevating Legrest



Standard For Models

KM-5000
KM-CP33

Optional For Models

S-Ergo 125
S-Ergo 115
KM-BT10
KM-7520
Econ 800/Econ 805

Axle traveling Elevating Legrest



Standard For Models

MVP 502

Optional For Models

Flexx S
Flexx
Flexx Tall
Flexx HD
VIP 515

Amputee Footrest



Standard For Models

None

Optional For Models

S-Ergo 125
S-Ergo 115
Flexx S
Flexx
Flexx Tall
Flexx HD
VIP 515
MVP 502
KM-5000
Econ 800/Econ 805

Vertical Footrest



Standard For Models

KM-7520

Optional For Models

Flexx S
Flexx
VIP 515
MVP 502
KM-5000
KM-CP33

Anti-tippers



Standard For Models

S-Ergo 106
VIP 515
MVP 502
KM-5000
KM-7520
KM-CP33

Optional For Models

S-Ergo 125
S-Ergo 115
S-Ergo 105
Ergo Lite
Ergo Lite 2
Ergo Live
KM-AT60
KM-1500
KM-2500/KM-2500L
Econ 800/Econ 805

Swing-away Anti-tipper



Standard For Models

KM-BT10

Optional For Models

Flexx S
Flexx
Flexx Tall
Flexx HD

Attendant Brake



Standard For Models

S-Ergo 125
S-Ergo 115
S-Ergo 105
S-Ergo 106
Ergo Lite
Ergo Lite 2
VIP 515
MVP 502
KM-BT10
KM-CP33
KM-1500

Optional For Models

Flexx S
Flexx
Flexx Tall
Flexx HD
KM-5000
KM-2500/KM-2500L
SOMA 100/SOMA 105
SM-100.3/SM-150.3
Econ 800/Econ 805

Tie-down Bracket



Standard For Models

None

Optional For Models

Flexx S
Flexx
Flexx Tall
Flexx HD
VIP 515

Pelvic Belt



Standard For Models

VIP 515
MVP 502
KM-BT10
KM-AA20

Optional For Models

S-Ergo 125
S-Ergo 115
S-Ergo 105
S-Ergo 106
Ergo Lite
Ergo Lite 2
Ergo Live
Flexx S
Flexx
Flexx Tall
Flexx HD
KM-5000
KM-7520
KM-CP33
KM-1500
KM-2500/KM-2500L
SOMA 100/SOMA 105
SM-100.3/SM-150.3
Econ 800/Econ 805

Harness



Standard For Models

None

Optional For Models

S-Ergo 125
S-Ergo 115
S-Ergo 105
S-Ergo 106
Ergo Lite
Ergo Lite 2
Flexx S
Flexx
Flexx Tall
VIP 515
MVP 502
KM-5000
KM-CP33
KM-1500
KM-2500/KM-2500L
SOMA 100/SOMA 105
SM-100.3/SM-150.3
Econ 800/Econ 805

Rigidifying Push Bar



Standard For Models

S-Ergo 106
Flexx HD
MVP 502
KM-BT10

Optional For Models

S-Ergo 125
S-Ergo 115
S-Ergo 105
Flexx S
Flexx
Flexx Tall
VIP 515

Foldable Headrest



Standard For Models

None

Optional For Models

S-Ergo 125
S-Ergo 115
S-Ergo 105
S-Ergo 106
Ergo Lite
Ergo Lite 2
Flexx S
Flexx
Flexx Tall
KM-1500
KM-2500/KM-2500L
Econ 800/Econ 805

Gel Cushion



For seat width : 16" x 16" / 18" x 16"

Standard For Models

None

Optional For Models


Flexx S
Flexx
VIP 515
KM-5000
KM-1500
KM-2500/KM-2500L
SOMA 100/SOMA 105
SM-100.3/SM-150.3
Econ 800/Econ 805

Lap Tray



Standard For Models

None

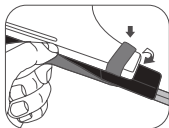

Optional For Models

S-Ergo 125	KM-5000
S-Ergo 115	KM-7520
S-Ergo 105	KM-CP33
S-Ergo 106	KM-1500
Flexx S	KM-2500/KM-2500L
Flexx	SOMA 100/SOMA 105
Flexx Tall	SM-100.3/SM-150.3
VIP 515	Econ 800/Econ 805
MVP 502	

Travel Bag



Standard For Models

None

Optional For Models

S-Ergo 125
S-Ergo 115
Ergo Lite
Ergo Lite 2
KM-2500/KM-2500L

Backrest Bag



For seat width : 16"-18"

Standard For Models

None

Optional For Models

S-Ergo 125	VIP 515
S-Ergo 115	KM-5000
S-Ergo 105	KM-1500
Ergo Lite	KM-2500/KM-2500L
Ergo Lite 2	SOMA 100/SOMA 105
Flexx S	SM-100.3/SM-150.3
Flexx	Econ 800/Econ 805
Flexx Tall	

Net Bag



For seat width : 16"-18"

Standard For Models

None

Optional For Models

S-Ergo 125	VIP 515
S-Ergo 115	MVP 502
S-Ergo 105	KM-5000
S-Ergo 106	KM-1500
Ergo Lite	KM-2500/KM-2500L
Ergo Lite 2	SOMA 100/SOMA 105
Flexx S	SM-100.3/SM-150.3
Flexx	Econ 800/Econ 805
Flexx Tall	

Cane Holder



Standard For Models

None

Optional For Models

S-Ergo 125	MVP 502
S-Ergo 115	KM-5000
S-Ergo 105	KM-7520
S-Ergo 106	KM-1500
Ergo Live	KM-2500/KM-2500L
Flexx S	SOMA 100/SOMA 105
Flexx	SM-100.3/SM-150.3
Flexx Tall	Econ 800/Econ 805
Flexx HD	
VIP 515	

I.V. Pole Holder



Standard For Models

None

Optional For Models

S-Ergo 125	KM-5000
S-Ergo 115	KM-BT10
S-Ergo 105	KM-7520
S-Ergo 106	KM-CP33
Flexx S	KM-1500
Flexx	KM-2500/KM-2500L
Flexx Tall	SOMA 100/SOMA 105
Flexx HD	SM-100.3/SM-150.3
VIP 515	Econ 800/Econ 805
MVP 502	

Oxygen Tank Holder



Standard For Models

None

Optional For Models

Flexx S	KM-7520
Flexx	KM-CP33
Flexx Tall	KM-1500
Flexx HD	Econ 800/Econ 805
VIP 515	
MVP 502	
KM-5000	
KM-BT10	

Canopy



Standard For Models

None

Optional For Models

S-Ergo 125
S-Ergo 115
S-Ergo 105
S-Ergo 106

Cup Holder



Standard For Models

None

Optional For Models

S-Ergo 125	KM-2500/KM-2500L
S-Ergo 115	SOMA 100/SOMA 105
S-Ergo 105	SM-100.3/SM-150.3
S-Ergo 106	
Ergo Lite	
Ergo Lite 2	
Flexx S	
VIP 515	
KM-1500	

Series	Prodcut Name	Seat Width (cm)	Seat Depth (cm)	Caster Size (inch)	Rear Wheel Size (inch/QR)	Seat Height (cm)	Backrest Height (cm)	Armrest Height (cm)	Calf Length (cm)
S-Ergo Series	S-Ergo 125	40	43	7	14F	50	41	21	31-40
					24F				
					24Q				
		45	43	7	14F	50	41	21	31-40
					24F				
					24Q				
		50	43	7	14F	50	41	21	31-40
					24F				
					24Q				
	S-Ergo 115	40	43	7	14F	50	41	21	31-40
					24F				
					24Q				
		45	43	7	14F	50	41	21	31-40
					24F				
					24Q				
		50	43	7	14F	50	41	21	31-40
					24F				
					24Q				
	S-Ergo 105	40	43	7	14F	50	41	21	31-40
					24F				
					24Q				
		45	43	7	14F	50	41	21	31-40
					24F				
					24Q				
	S-Ergo 106	45	42	7	24F	50	43	23	31-40
					24Q		42		
Ergo Lite Series	Ergo Lite	40 45	42	6	14F	48	42	21	35-40
	Ergo Lite 2	40	42	6	14F	48	42	21	35-40
					20F				
					20Q				
		45	42	6	14F	48	42	21	35-40
					20F				
20Q									
Ergo Semi Active	Ergo Live	36	38	5/6	24Q	45-51	33-48	23	36-44
			41						
		41	41	5/6	24Q	45-51	33-48	23	36-44
			46						
		46	41	5/6	24Q	45-51	33-48	23	36-44
			46						

Backrest Angle (X°)	Seat Angle (X°)	Legrest Angle (X°)	Max. User Weight (kg)	Weight w/o Accessories (kg)	Total Weight (kg)	Overall L x W x H (cm)	Overall Collapsed L x W x H (cm)
-	-	-	115	11	12.9	91*57*91	70*30*74
				12.4	14.3	99*58*91	80*31*74
				8.9		99*58*91	80*31*74
-	-	-	115	11.2	13.1	91*62*91	70*30*74
				12.6	14.5	99*63*91	80*31*74
				9.1		99*63*91	80*31*74
-	-	-	115	11.4	13.3	91*67*91	70*30*74
				12.8	14.7	99*68*91	80*31*74
				9.3		99*68*91	80*31*74
-	-	-	115	10.3	12.1	91*57*91	70*30*74
				11.7	13.5	99*59*91	80*31*74
				8.2		99*59*91	80*31*74
-	-	-	115	10.6	12.4	91*62*91	70*30*74
				12	13.8	99*63*91	80*31*74
				8.5		99*63*91	80*31*74
-	-	-	115	10.9	12.7	91*67*91	70*30*74
				12.3	14.1	99*68*91	80*31*74
				8.8		99*68*91	80*31*74
-	-	-	115	11.8	11.8	99*57*91	88*29*74
				13.2	13.2	102*58*91	98*31*74
				9.7		102*58*91	98*31*74
-	-	-	115	12.2	12.2	99*62*91	88*29*74
				13.6	13.6	102*63*91	98*31*74
				10.1		102*63*91	98*31*74
5,10,15	-	-	115	14.5	14.5	104*63*90	98*31*74
				11			
-	-	-	100	8.5	8.5	94*57*91	86*27*75
				8.7	8.7	94*62*91	
-	-	-	100	7.9	9.2	93*57*91	68*27*75
				9.2	10.5	93*60*91	78*31*75
				5.8		93*60*91	78*31*75
-	-	-	100	8.1	9.4	93*62*91	68*27*75
				9.4	10.7	93*65*91	78*31*75
				6		93*65*91	78*31*75
(-5),0,5	-	-	125	8.7	11.7	Min. 89*51*83 Max. 89*51*97 Min. 92*51*83 Max. 92*51*97	Min. 89*32*83 Max. 89*32*97
(-5),0,5	-	-	125	8.9	11.9	Min. 89*56*83 Max. 89*56*97 Min. 94*56*83 Max. 94*56*97	Min. 89*32*83 Max. 89*32*97
(-5),0,5	-	-	125	9.1	12.1	Min. 89*61*83 Max. 89*61*97 Min. 94*61*83 Max. 94*61*97	Min. 89*32*83 Max. 89*32*97

Series	Prodcut Name	Seat Width (cm)	Seat Depth (cm)	Caster Size (inch)	Rear Wheel Size (inch/QR)	Seat Height (cm)	Backrest Height (cm)	Armrest Height (cm)	Calf Length (cm)
Flexx Series	Flexx	38/40.5/43/45.5/48/50.5	41/43.5/46	5/6/7/8	14/20/22/24	39-53	43/51/38-53	21-29	38-48
	Flexx Tall	40.5/43/45.5/48	48.5/51	6/7/8	14/22/24	45.5-53	43/51/38-46	21-29	38-48
	Flexx HD	51/56/61	41/46/51	6/7/8	14/22/24	45.5-53	43/51/38-46	21-29	38-48
	Flexx S	38/40.5/43/45.5/48/50.5	41/43.5/46	5/6/7/8	14/20/22/24	39-53	43/51/38-53	22-24	38-48
Tilt-in-Space Series/Recline	VIP 515	40	42	7	14F	46	49	21-29	41-46
					20F				
					20Q				
		45	42	7	14F	46	49	21-29	41-46
					20F				
					20Q				
	MVP 502	37	42	7	14F	51	51	23-29	42-49
					22F				
		40	42	7	14F	51	51	23-29	42-49
					22F				
					22Q				
		42	42	7	14F	51	51	23-29	42-49
					22F				
		45	42	7	14F	51	51	23-29	42-49
					22F				
					22Q				
		45	46	7	14F	51	51	23-29	42-49
					22F				
					22Q				
		50	46	7	14F	51	51	23-29	42-49
					22F				
					22Q				
	KM-5000	41	43	8	14F	51	45	22	41-48
					24F				
					24Q				
		46	43	8	14F	51	45	22	41-48
					24F				
					24Q				
		46	45	8	14F	51	45	22	41-48
					24F				
					24Q				
		51	45	8	14F	51	45	22	41-48
					24F				
					24Q				
		51	45	8	14F	51	45	22	41-48
					24F				
					24Q				

Backrest Angle (X°)	Seat Angle (X°)	Legrest Angle (X°)	Max. User Weight (kg)	Weight w/o Accessories (kg)	Total Weight (kg)	Overall L x W x H (cm)	Overall Collapsed L x W x H (cm)
(-4)-16	3-7	-	130	Min. 9.3	Min. 14.2 (14F) Min. 14.7 (24Q)	Min. 102*56*76 Max. 107*69*106	Collapsed Width 34
(-4)-16	3-7	-	130	Min. 10.6	Min. 15.2 (14F) Min. 15.7 (24Q)	Min. 115*59*84 Max. 118*67*104	Collapsed Width 34
(-4)-16	3-7	-	170	Min. 10.6	Min. 16.5 (14F) Min. 17.0 (24Q)	Min. 108*71*84 Max. 118*81*104	Collapsed Width 35
(-4)-16	3-7	-	130	Min. 9.3	Min. 14.2 (14F) Min. 14.7 (24Q)	Min. 102*56*76 Max. 107*69*106	Collapsed Width 34
-	0-35	-	115	15.7	18.9	108*62*99(134)	73*34*76
				15.9	19.1	108*67*99(134)	78*38*76
				12.9			
-	0-35	-	115	16.6	19.8	108*66*99(134)	73*34*76
				16.7	19.9	108*71*99(134)	78*38*76
				13.7			
0-67	-	0-70	115	14.9	19.9	118*54*96(134)	85*30*76
				15.4	20.5	118*59*96(134)	85*36*76
0-67	-	0-70	115	15	20	118*57*96(134)	85*30*76
				15.9	21	118*61*96(134)	85*36*76
				12.9			
0-67	-	0-70	115	15.1	20.1	118*59*96(134)	85*30*76
				16.4	21.5	118*63*96(134)	85*36*76
0-67	-	0-70	115	15.2	20.2	118*62*96(134)	85*30*76
				16.8	21.9	118*65*96(134)	85*36*76
				13.1			
0-67	-	0-70	115	15.3	20.3	118*62*96(134)	85*30*76
				16.9	22	118*65*96(134)	85*36*76
				13.2			
0-67	-	0-70	115	15.8	20.8	118*67*96(134)	85*30*76
				17.6	22.7	118*70*96(134)	85*36*76
				13.7			
5-75	-	0-75	115	15	20.6	127*58*93(129)	91*34*93
				16.3	21.9	127*62*93(129)	91*39*93
				13.3			
5-75	-	0-75	115	15.4	21	127*63*93(129)	91*34*93
				16.7	22.3	127*67*93(129)	91*39*93
				13.7			
5-75	-	0-75	115	15.5	21.1	127*63*93(129)	91*34*93
				16.8	22.4	127*67*93(129)	91*39*93
				13.8			
5-75	-	0-75	115	15.8	21.4	127*68*93(129)	91*34*93
				17.1	22.7	127*72*93(129)	91*39*93
				14.1			
5-75	-	0-75	115	15.9	21.5	127*68*93(129)	91*34*93
				17.2	22.8	127*72*93(129)	91*39*93
				14.2			

Series	Prodcut Name	Seat Width (cm)	Seat Depth (cm)	Caster Size (inch)	Rear Wheel Size (inch/QR)	Seat Height (cm)	Backrest Height (cm)	Armrest Height (cm)	Calf Length (cm)
Heavy Duty	KM-BT10	56	46/51	7/8	24Q	39-50	41-51	25	32-37
		61							
		66							
		66	56/61	7/8	24Q	39-50	41-51	25	32-37
		71							
		76							
Airport	KM-AA20	39	42	6	24Q	49	42	19	32
Pediatric Series	KM-7520	28	28	6	12F	46	41	17-23	29
		30.5	30.5						
		34	34						
		28	28	6	20Q	46	41	17-23	29
		30.5	30.5						
		34	34						
	KM-7520 Recline	28	28	6	20Q	46	41	17-23	29
		30.5	30.5						
		34	34						
	KM-CP33	20.5-38 (Interval 2.5)	23-40.5 (Interval 2.5)	6	14F	45	25-56	15-33	38-46
					20F				
Active Series	KM-AT60	38-45.5 (Interval 2.5)	35.5-50.5 (Interval 2.5)	5	24Q	46-58	31-48	-	33-48
Standard Series	KM-1500	39	42	8	24F	48	38	21	32-43
		45							
	KM-2500	41	40	6	14F	47	39	21	32-43
		47							
	KM-2500L	41	40	6	20F	47	39	21	32-43
					20Q				
		47	40	6	20F	47	39	21	32-43
					20Q				
	KM-TV20.2	41	47	6	8F	53	47	28	-
		46							

Backrest Angle (X°)	Seat Angle (X°)	Legrest Angle (X°)	Max. User Weight (kg)	Weight w/o Accessories (kg)	Total Weight (kg)	Overall L x W x H (cm)	Overall Collapsed L x W x H (cm)
(-5)-15	-	-	250	22.1	30.2	Min. 118*78*80 Max. 123*88*101	Min. 90*44*80 Max. 95*44*101
				22.9	31		
				23.7	31.8		
(-5)-15	-	-	250	24.5	32.6	Min. 123*88*80 Max. 128*98*101	Min. 95*44*80 Max. 100*44*101
				25.4	33.5		
				26.2	34.3		
-	-	-	100	11.5	16.2	99*58*91	99*58*91
-	-	-	100	12	14	76*49*91	57*31*84
				12.3	14.3	79*50*91	60*31*84
				12.6	14.6	82*51*91	63*31*84
-	-	-	100	12	14.5	81*51*91	62*33*84
				12.3	14.8	84*52*91	65*33*84
				12.6	15.1	87*53*91	68*33*84
0-60	-	0-75	100	13.5	17	85*51*91	62*33*84
				13.8	17.3	88*52*91	65*33*84
				14.1	17.6	91*53*91	68*33*84
0-30	0-40	0-75	80	Min. 9	Min. 11	Min. 83*43*65 Max. 106*61*96	Collapsed Width 40
				Min. 9.5	Min. 11.5		Collapsed Width 38
-	0,3,5	-	100	Min. 8	Min. 10.5	Min. 88*62*77 Max. 103*70*106	-
-	-	-	100	13.5	13.5	100*63*87	91*31*72
				13.7	13.7	100*68*87	
-	-	-	100	10.9	10.9	89*53*86	72*25*68
				11.1	11.1	89*60*86	
-	-	-	100	11.7	11.7	89*56*86	82*28*68
				8.8			
-	-	-	100	11.9	11.9	89*63*86	82*28*68
				9			
-	-	-	100	8.7	8.7	91*55*92	61*35*71
				8.9	8.9	91*60*92	

Series	Prodcut Name	Seat Width (cm)	Seat Depth (cm)	Caster Size (inch)	Rear Wheel Size (inch/QR)	Seat Height (cm)	Backrest Height (cm)	Armrest Height (cm)	Calf Length (cm)
Entry Level Series	SOMA 100	40	42	7	16F	50	38	22	41-45
					22F				
		45	42	7	16F	50	38	22	41-45
					22F				
	SOMA 105	40	42	7	16F	50	38	22	41-45
					22F				
		45	42	7	16F	50	38	22	41-45
					22F				
	SM-100.3	43	42	7	16F	51	38	23	41-44
					22F				
	SM-150.3	43	42	7	16F	51	38	23	41-44
					22F				
	ECON800	41	43	8	16F	48	43	23	39-48
					24F				
					24Q				
		46	43	8	16F	48	43	23	39-48
					24F				
					24Q				
		51	43	8	16F	48	43	23	39-48
					24F				
					24Q				
	ECON805	41	43	8	16F	48	43	23	39-48
					24F				
					24Q				
		46	43	8	16F	48	43	23	39-48
					24F				
					24Q				
		51	43	8	16F	48	43	23	39-48
					24F				
					24Q				

Backrest Angle (X°)	Seat Angle (X°)	Legrest Angle (X°)	Max. User Weight (kg)	Weight w/o Accessories (kg)	Total Weight (kg)	Overall L x W x H (cm)	Overall Collapsed L x W x H (cm)
-	-	-	100	10.4	10.4	93*65*87	93*28*87
				11.1	11.1	97*65*87	97*28*87
-	-	-	100	10.6	10.6	93*68*87	93*28*87
				11.3	11.3	97*68*87	97*28*87
-	-	-	100	10.8	10.8	93*65*87	93*28*73
				11.5	11.5	97*65*87	97*28*73
-	-	-	100	11	11	93*68*87	93*28*73
				11.7	11.7	97*68*87	97*28*73
-	-	-	100	10.6	10.6	95*63*86	95*25*86
				11.3	11.3	99*64*88	99*29*88
-	-	-	100	10.8	10.8	95*63*86	95*25*73
				11.5	11.5	99*64*88	99*29*73
-	-	-	115	10.5	13	101*55*94	77*28*94
				12.1	14.6	111*59*94	87*32*94
				9.1			
-	-	-	115	10.7	13	101*60*94	77*28*94
				12.3	14.6	111*64*94	87*32*94
				9.3			
-	-	-	115	10.5	13.2	101*55*94	77*28*94
				12.1	14.8	111*59*94	87*32*94
				9.1			
-	-	-	115	10.9	13.4	101*65*94	77*28*75
				12.5	15	111*69*94	87*32*75
				9.5			
-	-	-	115	11	13.5	101*60*94	77*28*75
				12.6	15.1	111*64*94	87*32*75
				9.6			
-	-	-	115	11.2	13.7	101*65*94	77*28*75
				12.8	15.3	111*69*94	87*32*75
				9.8			

karma

A Better Wheelchair, A Better Fit

KARMA MEDICAL PRODUCTS CO., LTD.

No.2363,Sec.2,University Rd.,Min-Hsiung Shiang,Chia-Yi 621,Taiwan

T +886-5-206 6688 F +886-5-206 7788

Email: globalsales@karma.com.tw

www.KarmaMedical.com



ISO 9001

Power Wheelchairs

C O L L E C T I O N



karma

A Better Wheelchair, A Better Fit



Karma

"I can only truly feel alive while my wheels are spinning..."

These words from a wheelchair user express the importance of mobility and that's the importance we are committed to support. Such brave people deserve a wheelchair that suits their needs.

Entering the third decade of operation, "Fit the Needs of the Users" is always the main belief of Karma Medical Products Co., Ltd. The three major concepts of our innovation are: product development with specific features to correspond to users' medical needs, cutting edge ergonomic design and environmental-friendliness. We provide the very best wheelchairs and help each customer start a new life with greater freedom and comfort.

In the past 25 years, we have built up over 4,000 dealer-services successfully in more than 40 countries throughout 5 continents, supporting millions of wheelchair users and their families, enriching the lives of physically challenged people and sharing joyful and touching moments all around the world.

Global Marketing and Service Network



EUROPE

- U.K.
- Spain
- Belarus
- Benelux
- Czech
- Denmark
- Finland
- France
- Germany
- Greece
- Iceland
- Norway
- Poland
- Russia
- Spain
- Slovakia

ASIA

- Beijing
- Shanghai
- Huzhou
- Taiwan
- Thailand
- India
- Hong Kong
- Indonesia
- Japan
- Korea
- Malaysia
- Philippines
- Singapore
- Thailand
- Vietnam

AMERICA

- Argentina
- Canada
- Colombia
- Chile
- Mexico
- U.S.

OCEANIA

- Australia
- New Zealand

WEST ASIA

- Bahrain
- Iran
- Israel
- Kuwait
- Oman
- Qatar
- Saudi Arabia
- Turkey
- UAE

AFRICA

- Egypt
- South Africa



karma

INDEX

OUTDOOR/INDOOR

02
EVO ALTUS

04
MID LECTUS

06
EVO LECTUS

10
LEON

OUTDOOR

12
SABER

INDOOR/OUTDOOR

14
Standing Wheelchair
Ergo Stand / SME

18
eFlexx

20
Morgan

22
Blazer

24
KP-25.2

26
Ergo Nimble/
KP-10.3 CPT

28
SP-100

Seating System

30
Seating System for
Morgan/Leon

32
KISS Rehab Seat

Accessories

34

Specification

36

PRODUCTS MAP



EVO ALTUS

Stand Up, Stand Out

The new high-end standing wheelchair gives you freedom and independence in mobility. The standing and multiple power functions combined with a very low seat height of 41.5 cm and up to 80Ah batteries make EVO Altus the most unique standing wheelchair in the world!



Chest Support
An optional accessory that can protect the user when standing.



Low Seat Height
The 41.5 cm seat height makes it possible to sit at almost any table and it is also easy to travel in a small van or taxi.

Support Wheels
The powered support wheels move down automatically when standing function is activated.

Relax Mode
Power recline and power elevating legrest functions allow the user to lie down straight as total flat position.

45° Power Tilt-in-space
Releases pressures at buttocks and stabilizes body position while driving downhill.

Power Elevating
Up to 40cm lifting function makes users reach things in high place easily. In daily social activities, users can have eye contact with others.



Standing Functions

Three standing modes controlled by an on-board computer



1. Sit-to-stand
Provides the most natural way to stand up. Also, the Altus seat helps you stand up in a safe and stable way.



2. Lay-to-stand
Offers a good preparation before standing up, and serves as a fine training tool to stretch out sedentary tendons without burdening your body.



3. 45° lay-to-stand
For the users who feel from lying to standing is too unnatural, this mode is a good compromise indeed.

OPTIONS & ACCESSORIES

- Chest Support
- Transfer Handles
- Push Handles
- Rear Bumper Roller
- Parallel Swing-Away-Down Controller Bracket
- Table Width Flip Around Control
- 3D Adjustable Headrest
- 3D Lateral Support
- 3D Hip Support
- 3D Calf Support
- Positioning Belt
- Positioning Belt Retractable
- LED Lighting System Include Indicators

Special controls & individual adoptions are available on request

SPECIFICATION

Seat Width (cm)	42 - 52
Front/Rear Wheel (inch)	14 / 9
Power Functions	Standing Tilt 0°-45° Bio-mechanical Recline (-5°)-90° Bio-mechanical Legrest: (-5°)-90° Lift: 0-40 cm
Controller	Rnet
Motor	320W
Battery	50Ah*2 / 80Ah*2
Min. Turning Radius (cm)	75
Min. Ground Clearance (cm)	7.5
Max. Safe Slope	10°
Max. Obstacle Height (cm)	7
Max. Speed (km/h)	6,10,12
Max. Driving Range (km)	25 / 45
Max. User Weight (kg)	140
w/ Batteries / w/o Batteries (kg)	215(80Ah)/ 165
Overall Length x Width x Heigh (cm)	120*61*112

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 36-37.
-Battery is not part of the product, might change as per local country specification and availability.

Create ultimate experience from the heart. Through the distinguished mechanical design, we make this mid-wheel-drive wheelchair with the best maneuverability and high performance. Furthermore, its full-scale adjustable seating system with power position changes will provide you the most comfortable experience than ever. Perform a wonderful spin with MID Lectus!



LECTUS Seat


Release body pressure with the Lectus positioning seating system, its multiple adjustments and smooth texture support a pleasurable sitting.



MID Chassis

With silent chassis and 6 wheel linkage suspension system, MID Lectus brings you the utmost comfort, stable and smooth maneuverability for both indoor and outdoor use.

6 WHEEL SUSPENSION SYSTEM

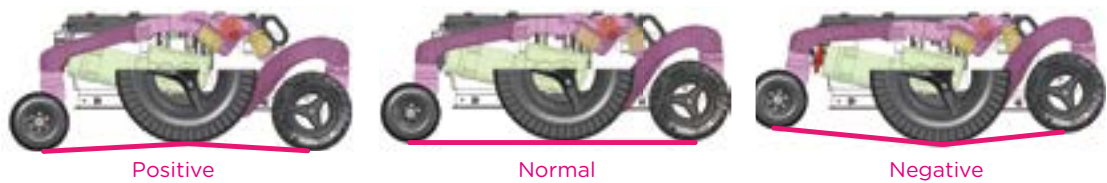

Designed with linked suspension on the front and rear caster, it can evenly distribute the driving impact and enable to overcome curbs up to 7.5 cm.



45° Power Tilt-in-space

Release the pressure at buttocks and stabilize the body position while driving downhill.



30cm Power Elevating System

Easy to reach higher spaces.

180° Power Recline + 180° Power Elevating Legrest

Total flat position allows the user to lie down straight.



Advanced Driving Range

Larger battery box that can carry 80Ah battery, with embedded battery design provides tool-free assembly for convenient maintenance. Go farther from now on!

OPTIONS & ACCESSORIES

- Transfer Handles
- Push Handles
- Rear Bumper Roller
- Parallel Swing-Away-Down Controller Bracket
- Table Width Flip Around Control
- 3D Adjustable Headrest
- 3D Lateral Support
- 3D Hip Support
- 3D Calf Support
- Positioning Belt
- Positioning Belt Retractable
- LED Lighting System Include Indicators

Special controls & individual adoptions are available on request

SPECIFICATION

Seat Width (cm)	42 - 52
Front/Middle/Rear Wheel (inch)	8 / 14 / 6
	Tilt 0°-45°
Power Functions	Bio-mechanical Recline (-5°)-90° Bio-mechanical Legrest: (-5°)-90° Lift: 0-30cm
Controller	Rnet
Motor	350W
Battery	50Ah*2 / 80Ah*2
Min. Turning Radius (cm)	62
Min. Ground Clearance (cm)	9
Max. Safe Slope	10°
Max. Obstacle Height (cm)	7.5
Max. Speed (km/h)	6,10,12
Max. Driving Range (km)	24/36
Max. User Weight (kg)	140
w/ Batteries / w/o Batteries (kg)	190(80Ah)/ 133
Overall Length x Width x Heigh (cm)	110*61.5*104

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 36-37.
-Battery is not part of the product, might change as per local country specification and availability.




EVO LECTUS represents an evolution from years of experience in wheelchair assessments. It's born to suit the exact specification of each market, and can be customized to fit the needs of every single client.



LECTUS Seat

The full-scale adjustable system can be set up to fit the individual needs, including fully multi-adjustable seat, armrests and legrests.

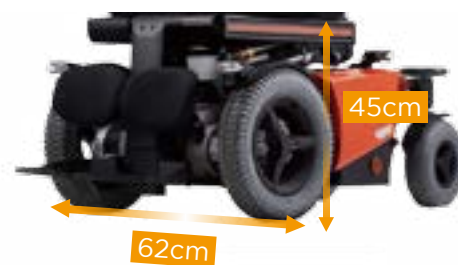

EVO Chassis

Equipped with a 4-wheel superb full suspension system which guarantees powerful and smooth driving performances in all terrains.

With 45cm low seat height, 62cm in total width, and 88cm in length after the footplates are flipped up, EVO Lectus provides an agile using experience.

40cm Power Elevating

Allows seat height from 45 cm up to 85cm.



180° Power Recline + 180° Power Legrest

Total flat position enables the user to lie down straight.



50° Power Tilting

Dynamic COG design allows tilting up to 50 degree.



Modular Design

1.Modular Power Settings
Including dynamic COG tilt, biomechanical recline and legrest elevating, and seat lift.

2.Multi-adjustable
Armrest (5 axis) & Footrest (3 axis)

3.Every function can be upgraded by module.

Advanced Electronics
EVO LECTUS is compatible with R-net to provide various possibilities.



OPTIONS & ACCESSORIES



- Transfer Handles
- Push Handles
- Rear Bumper Roller
- Parallel Swing-Away-Down Controller Bracket
- Table Width Flip Around Control
- 3D Adjustable Headrest
- 3D Lateral Support
- 3D Hip Support
- 3D Calf Support
- Positioning Belt
- Positioning Belt Retractable
- LED Lighting System Include Indicators

Special controls & individual adoptions are available on request

A Chair That
Evolves with You



SPECIFICATION

Seat Width (cm)	42 - 52
Front/Rear Wheel (inch)	14 / 9
Power Functions	Tilt 0°-50° Bio-mechanical Recline (-5°)-90° Bio-mechanical Legrest: (-5°)-90° Lift: 0-40cm
Controller	Rnet
Motor	320W
Battery	80Ah*2
Min. Turning Radius (cm)	75
Min. Ground Clearance (cm)	7
Max. Safe Slope	10°
Max. Obstacle Height (cm)	7
Max. Speed (km/h)	6/10
Max. Driving Range (km)	40
Max. User Weight (kg)	140
w/ Batteries / w/o Batteries (kg)	210/160
Overall Length x Width x Heigh (cm)	120*62*112

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 36-37.
-Battery is not part of the product, might change as per local country specification and availability.

LEON®



The Leon is designed to exploit the advantages of a rear-wheel-drive chair for outdoor use. It is geared for speed without compromising on stability or control. The exceptionally narrow wheelbase allows for agile maneuverability. The Leon is a superior power chair for even the most active users.



Superior Riding Comfort

Equipped with 4-wheel independent suspension for extra comfort and ultimate riding quality.



Excellent Stability

- High performance rear-wheel drive with continuous acceleration.
- 450W 4-pole powerful motors bring speed without compromising on stability or control.

Multiple Choices of Seats

- 17"/19" Sling Seat
 - 18"/20" Captain Seat
 - KISS Rehab Seat
- For more information, please refer to page 30-31.



Flip-back Armrest with Lock Mechanism

Easy to do lateral transfer, and better driving stability.



High Cruising Power

The optional 80Ah battery provides an extended driving range of up to 28 miles or 45km for long distance outdoor use.

*Exceed Your Expectation
Excel in Navigation*



OUTDOOR/INDOOR LEON

OPTIONS & ACCESSORIES

- LED Lighting System
- Swing Away Controller Bracket
- Contour Adjustable Headrest (for Sling Seat only)
- Elevating Legrest
- Hand bar

PRODUCT VIDEO LINK



SPECIFICATION

	Sling Seat	Captain Seat
Seat Width (cm)	43 / 48	46 / 51
Front/Rear Wheel (inch)	9 / 14	9 / 14
Power Functions	Tilt:5°-23° & Recline: 0-32°	-
Controller	VR2	VR2
Motor	450W	450W
Battery	50Ah*2/80Ah*2	50Ah*2/80Ah*2
Min. Turning Radius (cm)	78	78
Min. Ground Clearance (cm)	6.5	6.5
Max. Safe Slope	8°	8°
Max. Obstacle Height (cm)	5	5
Max. Speed (km/h)	6/10/12	6/10/12
Max. Driving Range (km)	30/45	30/45
Max. User Weight (kg)	140	140
w/ Batteries / w/o Batteries (kg)	108/78	100/70
Overall Length x Width x Height (cm)	106*62*130	106*62*123

±1cm ±0.5kg

-Karma reserves the right to modify information herein without further notice.

-For more details, please refer to the Specification Pages No. 36-37.

-Battery is not part of the product, might change as per local country specification and availability.



In the winter of 2004, a Korean man with his SABER (KP-45.3) traversed from Los Angeles to New York, total 17 states, 10,000 km. It took him 5 months to achieve this remarkable journey. SABER is known for its robustness and endurance, the latest SABER(KP-45.5) has maintained the outstanding performance and upgraded with ISR suspension. It is the best companion for an explorer.



Independent Suspension Response (ISR)

Allows the absorbers to accurately predict and avoid road shock, significantly improving stability and comfort.

Off-road Trial

Successfully travelled over 1,000 km of the road test, passed the 6,667 times of the drop test, demonstrating a superb off-road safety.



Splendid LED Lighting

LED lighting control system enhances driving safety at night.



Optional Power Tilt and Recline Function

Provide superior pressure relief and postural support.



Flip-back Armrest with Lock Mechanism

Make it easier to transfer into the chair.



The Best Companion for an Explorer

OUTDOOR SABER

OPTIONS & ACCESSORIES

- Swing Away Controller Bracket
- Elevating Legrest
(Standard Specification for SABER TR)
- Hand bar
- Red ISR system for users from 110kgs to 135kgs
(Standard Specification for SABER TR)



SPECIFICATION

	SABER	SABER TR
Seat Width (cm)	45 / 51	51
Front/Rear Wheel (inch)	9.5 / 16	9.5 / 16
Power functions	-	Tilt: 9°-47° / Recline: 0°-25°
Controller	VR2	VR2
Motor	400W	400W
Battery	50Ah*2/80Ah*2	50Ah*2
Min. Turning Radius (cm)	88	97
Min. Ground Clearance (cm)	10	10
Max. Safe Slope	10°	10°
Max. Obstacle Height (cm)	8	8
Max. Speed (km/h)	13	13
Max. Driving Range (km)	38km(50Ah)/57km(80Ah)	38km(50Ah)
Max. User Weight (kg)	135	135
w/ Batteries / w/o Batteries (kg)	102/73	113/82
Overall Length x Width x Height (cm)	117.5*65*120	120*65*126

±1cm ±0.5kg
 -Karma reserves the right to modify information herein without further notice.
 -For more details, please refer to the Specification Pages No. 36-37.
 -Battery is not part of the product, might change as per local country specification and availability.



Standing Series

A wheelchair is no longer just a mobility product; it is an important tool to join the crowd and return to the desirable social life.

Height can be presented a big problem for wheelchair users. Standing wheelchair is to enable wheelchair users to stand up with ease, to enjoy more independent living, and to acquire an equal opportunity for face-to-face interactions with friends and relatives.




Standing can activate the main circulation systems and help maintain physical and mental health.

- Activates circulation system and improves cardio-pulmonary function.
- Helps prevent decubitus sores.
- Improves bowel regularity and function.
- Improves bladder function and decreases urinary tract infections.
- Reduces contractures and ankylosis.
- Prevents osteoporosis.
- Improves quality and independence of life.



Power Recline Function
Provides a wide angle adjustment of (-10°) - 30° with excellent back support.

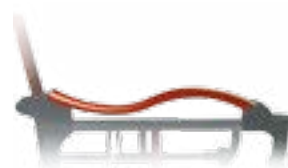

S-Ergo Seating System
The unique seating system can release pressure, stabilize posture, and reduce sliding.



Front-wheel Drive
Provides agility, easy movements, and steady standing experience. It is also easy to turn around in a small space.



Power Standing Wheelchair

OPTIONS & ACCESSORIES

- Lighting System
- Swing Away Controller Bracket
- Foldable & Adjustable Headrest
- H-shaped Belt
- Pelvic Belt
- Mud Guard
- Assistant Brake

SPECIFICATION

Seat Width (cm)	40 / 44 / 48
Front / Rear Wheel (inch)	12 / 6
Power Functions	Power stand & Recline:(-10°) - 30°
Controller	VR2
Motor	300W
Battery	36Ah*2/50Ah*2
Min. Turning Radius (cm)	67
Min. Ground Clearance (cm)	6.5
Max. Safe Slope	10
Max. Obstacle Height (cm)	5
Max. Speed (km/h)	9
Max. Driving Range (km)	28(50Ah)
Max. User Weight (kg)	135
w/ Batteries / w/o Batteries (kg)	89 / 57
Overall Length x Width x Height (cm)	102*63*98

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 36-37.
-Battery is not part of the product, might change as per local country specification and availability.



Manual Standing Wheelchair

OPTIONS & ACCESSORIES

- Foldable & Adjustable Headrest
- H-shaped Belt

SPECIFICATION

Seat Width (cm)	40 / 44 / 48
Seat Depth (cm)	44 - 50
Front/Rear Wheel (inch)	6 / 24
Seat Height (cm)	52
Backrest Height (cm)	42
Armrest Height (cm)	18-26
Legrest Length (cm)	36-42
Battery	2.9Ah*2
Max. User Weight (kg)	120
w/ Batteries / w/o Batteries (kg)	28.6 / 25.8
Overall Length x Width x Height (cm)	104*(59-67)*(96-143)

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 36-37.
-Battery is not part of the product, might change as per local country specification and availability.



Quick Release Rear Wheels
Reduce the weight and size for easy carrying and storage.




Four-point Touch Ground
Provides extra safety in standing position.

eFlexx



The eFlexx was created to provide the user with the best possible experience available, which is simple, light and efficient. Without the need of taking off the battery, eFlexx can be easily folded into a compact form for ease of transportation.

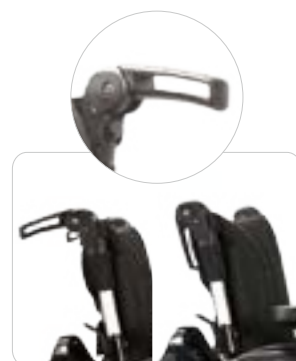

Foldable Push Handles

Allow minimal effort in folding and pushing the wheelchair. It also makes eFlexx look sleek and streamline.



Swing-away Controller System

Easier to approach the desk-top. The new safety feature prevents hands and fingers from pinching accidents.



in a folded state

Lightweight & Removable Battery

Battery installation can be easily done in just one simple step. Allows to charge anywhere and anytime.

- Standard 26.1Ah lithium battery provides an extended driving range.
- Optional 11.5Ah lithium battery makes travel friendly.



Benefit of Lithium battery !

Ecological friendly, highly efficient and long lasting. The battery weight (3.9 kg) is only 1/6 the weight of the traditional lead-acid batteries.



-4°/0°/4°/8°/12°/16°
Angle Adjustable Backrest

Easy to operate angle adjustment from -4° to 16°, giving users a more comfortable experience and flexible choices.



*Travel Freely
and Lightly.*

INDOOR/OUTDOOR eFlexx

OPTIONS & ACCESSORIES

- Boat Armpad
- Hemi Armpad
- Elevating Legrest
- Assistant Brake
(Standard for EFL-20)
- Tie-down Bracket
- Rigidifying Push Bar
(for 43,48cm seat width only)
- Cane Holder

SPECIFICATION

	EFL-14	EFL-20
Seat Width (cm)	38 / 43 / 48	40.5 / 43 / 45.5
Front/Rear Wheel (inch)	8 / 14	8 / 20
Controller	Linx	Linx
Motor	250W	250W
Battery	26.1 / 11.5	26.1 / 11.5
Min. Turning Radius (cm)	90	90
Min. Ground Clearance (cm)	7	7
Max. Safe Slope	6°	6°
Max. Obstacle Height (cm)	2.5	2.5
Max. Speed (km/h)	8	8
Max. Driving Range (km)	20(26.1Ah)	20(26.1Ah)
Max. User Weight (kg)	100	100
w/std Batteries / w/o Batteries (kg)	33 / 29	35.5 / 31.5
Overall Length x Width x Height (cm)	110*65*85	117.5*66*85

±1cm ±0.5kg

-Karma reserves the right to modify information herein without further notice.

-For more details, please refer to the Specification Pages No. 36-37.

-Battery is not part of the product, might change as per local country specification and availability.



Morgan is a mid-wheel-drive wheelchair with high performance to meet both indoor and outdoor requirement. Along with its modern and sleek design, Morgan is an excellent ride ideally for those who want a great all-round companion.



Excellent Indoor Maneuverability

With a tight turning radius of only 64 cm, the Morgan has great indoor control to drive in narrow spaces. It is easy to move around and navigate through doorways or into the lifts with a narrow width of just 60 cm.



Superior Intuitive Driving

Great intuitive driving is the main benefit carried by the Morgan. This is ideal for those who may have difficulty in driving powered devices or operating a joystick.

Multiple Choices of Seats

- 17"/19" Sling Seat
 - 18"/20" Captain Seat
 - KISS Rehab Seat
- For more information, please refer to page 30-31.

Flip-back Armrest with Lock Mechanism

Easy to do lateral transfer, and better driving stability.

Smooth, Stable and Safe

The six-wheel independent suspension system ensures a smooth and stable drive while riding on rough terrain.

A Great All-round Companion



OPTIONS & ACCESSORIES

- LED Lighting System
- Swing Away Controller Bracket
- Contour Adjustable Headrest (for Sling Seat only)
- Elevating Legrest
- Hand bar

PRODUCT VIDEO LINK



SPECIFICATION

	Sling Seat	Captain Seat
Seat Width (cm)	43 / 48	46 / 51
Front/Rear Wheel (inch)	8 / 14 / 6	8 / 14 / 6
Power Functions	Tilt:5°-23° & Recline: 0-32°	-
Controller	VR2	VR2
Motor	450W	450W
Battery	50Ah*2	50Ah*2
Min. Turning Radius (cm)	64	64
Min. Ground Clearance (cm)	8	8
Max. Safe Slope	8°	8°
Max. Obstacle Height (cm)	5	5
Max. Speed (km/h)	6/10	6/10
Max. Driving Range (km)	25	25
Max. User Weight (kg)	140	140
w/ Batteries / w/o Batteries (kg)	116/86	108/78
Overall Length x Width x Height (cm)	115*60*130	115*60*123

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 36-37.
-Battery is not part of the product, might change as per local country specification and availability.



The compact design of Blazer provides not only a hassle free outdoor experience, but also a great indoor agility. It offers essential configurations with premium adjustments that exceed its value.



Indoor Agility
With only 60cm and total width and 107cm total length, it is highly maneuverable in tight area.



Outdoor Performance
9" caster, 14" rear wheel with suspension and mudguard offer an excellent outdoor experience and ensure the safety under 8° slope and 6 cm curb environment.



A Compact and Powerful Solution

INDOOR/OUTDOOR Blazer

Modular Power Function Options and Mutiple Choices of Seats



Optional Power Tilt-in-space
Provides instant pressure relief.



17"/19" Sling Seat
· Premium seat with adjustable seat depth (41-46 cm) and backrest angle ((-1)° - 14°)
· Compatible with Power Seat Tilt system.



18"/20" Captain Seat
· Comfortable seat with manual recline angle adjustment 0°-65°

OPTIONS & ACCESSORIES

- LED Lighting System
- Swing Away Controller Bracket
- Contour Adjustable Headrest (for Sling Seat only)
- Elevating Legrest
- Assistant Brake
- Hand bar
- Tie-down bracket

PRODUCT VIDEO LINK



SPECIFICATION

	Sling Seat	Captain Seat
Seat Width (cm)	43 / 48	46 / 51
Front/Rear Wheel (inch)	9 / 14	9 / 14
Power Functions	Tilt 4°-25°	-
Controller	VR2	VR2
Motor	300W	300W
Battery	50Ah*2	50Ah*2
Min. Turning Radius (cm)	82	82
Min. Ground Clearance (cm)	7	7
Max. Safe Slope	8°	8°
Max. Obstacle Height (cm)	6	6
Max. Speed (km/h)	10	10
Max. Driving Range (km)	30	30
Max. User Weight (kg)	140	140
w/ Batteries / w/o Batteries (kg)	94/64	100/70
Overall Length x Width x Height (cm)	107*60*93	107*63*128

±1cm ±0.5kg

-Karma reserves the right to modify information herein without further notice.

-For more details, please refer to the Specification Pages No. 36-37.

-Battery is not part of the product, might change as per local country specification and availability.

KP-25.2

KP-25.2 is suitable for mid-short distance journey. Featured with light weight and foldable frame, compact size and easily removal battery pack, KP-25.2 is portable and convenient for your daily use.



Highly Configurable Meets High Adaptability

Compact medium-size wheelchair with overall width of 59 cm.



Height Adjustable Armrest

The armrest height is with four-level adjustment which can be adjusted according to users' body requirements. It allows users to relax shoulder muscles and offers better support for elbow.



Tool-free Quick Folding

Light weight, foldable frame and quick release battery removal design make transport even more easier.



Lateral Movement

It is convenient for lateral movement with flip-up armrest and swing-away detachable footrests.

Starting a New Life
with a New Scope



OPTIONS & ACCESSORIES

- Swing Away Controller Bracket
- Quick Release Attendant Controller
- Foldable & Adjustable Headrest
- Elevating Legrest
- Hand bar
- Pneumatic tires for rear wheel

SPECIFICATION

Seat Width (cm)	36 / 41 / 46 / 51
Front/Rear Wheel (inch)	8 / 12
Controller	VR2
Motor	200W
Battery	22Ah*2 / 36Ah*2 / 50Ah*2
Min. Turning Radius (cm)	82
Min. Ground Clearance (cm)	6
Max. Safe Slope	10°
Max. Obstacle Height (cm)	5
Max. Speed (km/h)	6
Max. Driving Range (km)	38(50Ah)
Max. User Weight (kg)	115 / 135 (w/ pne. tire)
w/ Batteries / w/o Batteries (kg)	70/35
Overall Length x Width x Heigh (cm)	110*(54-69)*96

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 36-37.
-Battery is not part of the product, might change as per local country specification and availability.

Ergo Nimble is an affordable rear wheel drive power wheelchair. It features the patented S-Ergo Seating System to promote a comfortable and good seating posture. Compact design makes it ideal for confined spaces mobility and short distance journey.



S-Ergo Seating System
The unique seating system can release pressure, stabilize posture, and reduce sliding.



Tool-free Angle Adjustment
Between working angle 5° and relaxing angle 19°.



Can be easily disassembled into 3 pieces
The tool-free disassemble system makes traveling easier. The maximum weight per piece is 26 kg, and they can easily be put in a car trunk.

OPTIONS & ACCESSORIES

- Swing Away Controller Bracket
- Quick Release Attendant Controller
- Assistant Brake

SPECIFICATION

Seat Width (cm)	44
Front / Rear Wheel (inch)	8 / 12
Controller	VR2
Motor	200W
Battery	22Ah*2
Min. Turning Radius (cm)	75
Min. Ground Clearance (cm)	9
Max. Safe Slope	8°
Max. Obstacle Height (cm)	3
Max. Speed (km/h)	6
Max. Driving Range (km)	12
Max. User Weight (kg)	120
w/ Batteries / w/o Batteries (kg)	50 / 35
Overall Length x Width x Heigh (cm)	103*62*94

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 36-37.

KP-10.3 CPT

Specifically designed to solve the problem of loading power wheelchairs into vehicles, our “KP-10.3 CPT” presents a revolutionary design that can be easily disassembled into three pieces without requiring any tools. Going out is no longer a challenge !



Comfortable Captain Seat
High quality seat offers extra riding comfort.



Height Adjustable and Detachable Armrest
Easy to fit and relax shoulder muscle without tools.



Carefree traveling
Disassembly and assembly are in simple steps without requiring any tools.

OPTIONS & ACCESSORIES

- Swing Away Controller Bracket
- Assistant Brake

SPECIFICATION

Seat Width (cm)	45
Front / Rear Wheel (inch)	8 / 12
Controller	VR2
Motor	200W
Battery	22Ah*2
Min. Turning Radius (cm)	75
Min. Ground Clearance (cm)	9
Max. Safe Slope	8°
Max. Obstacle Height (cm)	3
Max. Speed (km/h)	6
Max. Driving Range (km)	12
Max. User Weight (kg)	120
w/ Batteries / w/o Batteries (kg)	53.3/37.5
Overall Length x Width x Heigh (cm)	103*61*92

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 36-37.
-Battery is not part of the product, might change as per local country specification and availability.

SP-100

SP-100 is a highly cost effective lightweight power wheelchair, which is the best choice for the first-time user. Durable with chic design, SP-100 is the ultimate choice that can fulfill user’s need, and it is just beyond your expectations.



Easy Transfer
Flip-back armrests and foldable footplates, allows for easier transfer in and out of the chair.



Seat with Comfort and Easy to Maintain
Double-layered seat upholstery is removable and easy to clean.



Reinforced Double Crossbraces
Give you a comfortable and safety seating experience.



Lightweight Robust-structure Frame
The elliptical tubes create a modern and stylish look.



Best Choice for City Lohas
SP-100 is equipped with removable batteries. And, it can be folded into compact form even while the batteries are loaded. It fits the need of LOHAS that enjoy living freely.

Simplicity and Utility Design



INDOOR/OUTDOOR SP-100

OPTIONS & ACCESSORIES

- Direct-to-battery charging cable

SPECIFICATION

Seat Width (cm)	41/46
Front / Rear Wheel (inch)	8/16
Controller	Linx
Motor	100W
Battery	22Ah*2
Min. Turning Radius (cm)	81
Min. Ground Clearance (cm)	7
Max. Safe Slope	6°
Max. Obstacle Height (cm)	4
Max. Speed (km/h)	6
Max. Driving Range (km)	15
Max. User Weight (kg)	100
w/ Batteries / w/o Batteries (kg)	43 / 29
Overall Length x Width x Heigh (cm)	99*(56/61)*86

±1cm ±0.5kg
-Karma reserves the right to modify information herein without further notice.
-For more details, please refer to the Specification Pages No. 36-37.
-Battery is not part of the product, might change as per local country specification and availability.

Seating System

Adaptable
Products



1 Modular sling seating design improves comfort

The two power functions can be applied alone or in combination with each other depending on the users' preference and requirements.



Power Seat Tilt,
angle up to 23 degrees



Power Backrest Recline,
angle up to 32 degrees

2 Multiple choices of seats



17"/19" Sling Seat

Offer basic support

Other options are adaptable, such as air cushion, backrest, positioning assembly from aftermarket.



18"/20" Captain Seat

Shock attenuation

Comfortable and suitable for the user with better waist equilibrant.



KISS Rehab Seat

Expansible and adjustable

Suitable for progressive symptom.

- Premium seat with adjustable seat depth(41-46 cm) and backrest angle((-1)° - 11°).
- Compatible with **Power Seat Tilt** and **Power Backrest Recline** systems.

- Comfortable seat with manual recline angle adjustment
18" seat : 0°-65°
20" seat : 0°-45°

- Full range adjustment
- Compatible with **Power Tilt** with "Center of Gravity" weight shift, **Power Recline** with shearing force reduction design, **Power Elevating Legrest** and **Power Seat Lift**



Kiss

Designed to Conquer
Positioning Challenge


Karma Integrated Seating System offers numerous thoughtful functions that make daily life easier. The unique Push-2-Flip design is built in all supporting accessories to minimize the load of caregivers while transferring.

A fully adjustable seat frame, multiple power features and various accessories make KISS be the best choice in assisted mobility devices.

Fully Adjustable Seating System

KISS allows for width, depth, back height and angle adjustments to be made easily and quickly in the field, which can grow with the client.

	KISS Adult	KISS Junior
Seat Width	16"-20"	12"-16"
Seat Depth	16"-20"	12"-16"
Backrest Height	20"-22"	16"-18"

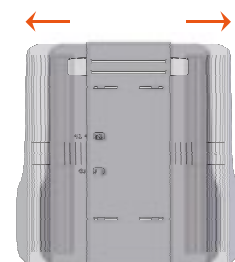

Backrest Pan



Seat Pan



Backrest Height



Seat Depth



Power Lift

Conquer the challenges of different height in daily life by up to 30cm lift.




Power Tilt

50° tilt provides perfect pressure distribution from bottom to back and enhances the posture stability. The advanced center of gravity system provides the best safety and stability.



Power Recline

With power elevating legrest and 170° shearing force reduction recline, KISS provides the best seating comfort and pressure relief.

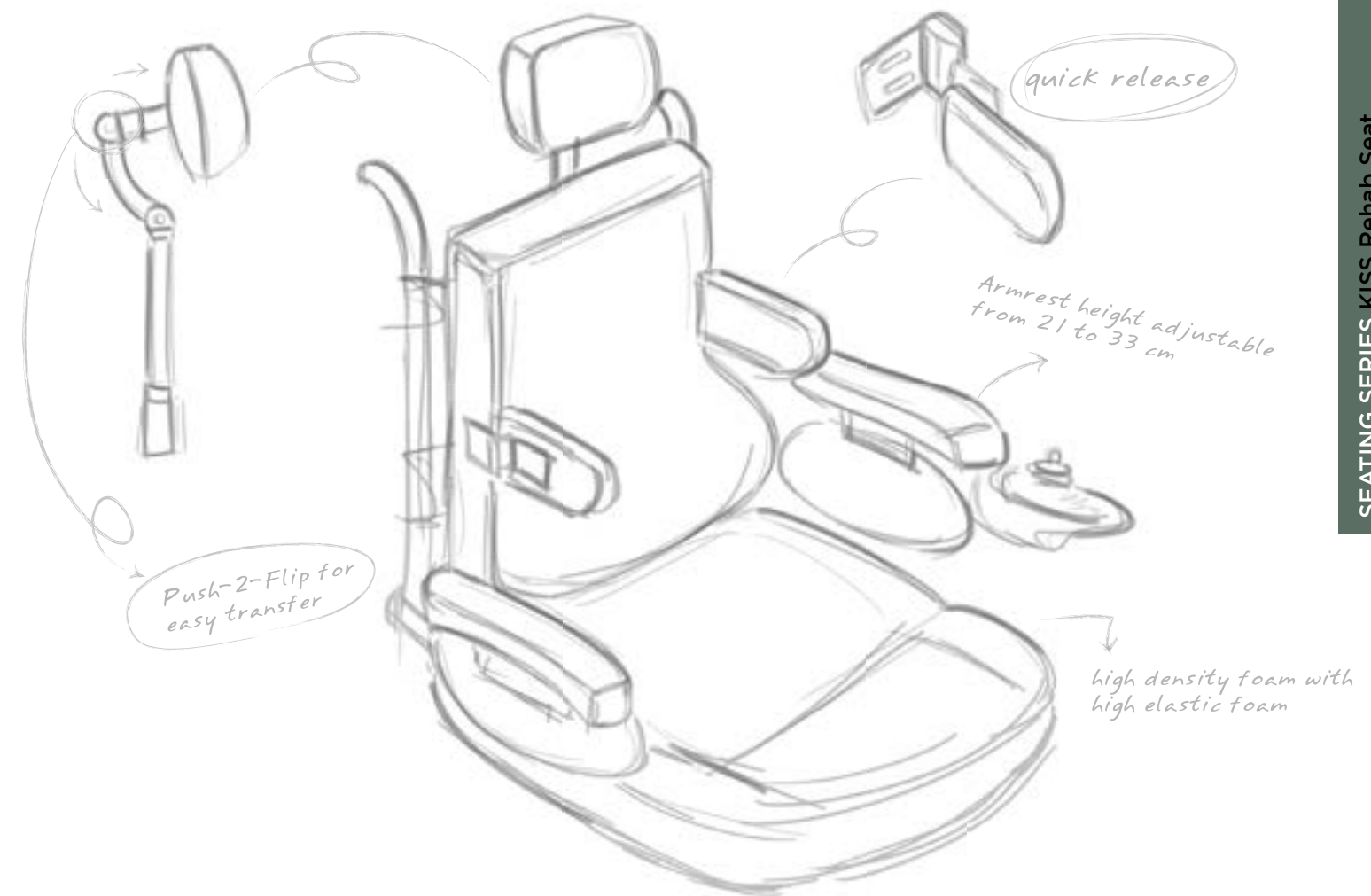


Modular Supporting Accessories

Modular supporting accessories with more than a hundred combinations for advanced positioning. Sliding rail for easy adaption and adjustment.



Trunk support

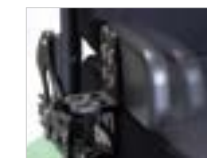

Hip support



Pommel

Push-2-Flip Design

With one simple button, caregivers can easily swing away supporting accessories and complete the transfer.



OPTIONS & ACCESSORIES

- Headrest
(Flip back, Curve adjustable)
- Swing-away Controller Bracket
- Armrest
(Flip-back, T-shape, Reclining)
- Armpad
(Ergo, Contour, Padded)
- Legrest
(Vertical, Standard, Elevating, Power elevating)
- Footplate
(Angle adjustable, Multi-adjustable)
- Modular Supporting Accessories
(Trunk, Hip, Pommel)

SPECIFICATION

	Kiss Adult	Kiss Junior
Seat Width (cm)	41 - 51	31 - 41
Seat Depth (cm)	41 - 51	31 - 41
Front Seat Height (cm)	44 - 49	44 - 49
Backrest Height (cm)	51 - 56	41 - 46
Seat Angle (X°)	0°, 3°, 6°/ -6°, -3°, 0°, 3°, 6°(T)	0°, 3°, 6°/ -6°, -3°, 0°, 3°, 6°(T)
Back Angle (X°)	0 - 30°	0 - 30°
Armrest Height (cm)	20 - 30	15 - 25
Legrest Length (cm)	32 - 46	28 - 37
Power Tilt (X°)	0 - 47	0 - 47
Power Recline (X°)	0 - 80	0 - 80
Power Legrest (X°)	0 - 80	0 - 80
Power Lift (cm)	30	30
Max. User Weight (kg)	140	100

±1cm ±0.5kg


-Karma reserves the right to modify information herein without further notice.

-For more details, please refer to the Specification Pages No. 36-37.

-Battery is not part of the product, might change as per local country specification and availability.

Accessories


Swing Away Controller Bracket



For Models

- LEON
- SABER
- SABER TR
- Ergo Stand
- MorgZan
- Blazer
- KP-25.2
- Ergo Nimble
- KP-10.3 CPT
- Kiss


Foldable & Adjustable Headrest



For Models

- Ergo Stand
- SME
- KP-25.2


Quick Release Attendant Controller



For Models

- KP-25.2
- Ergo Nimble


Contour Adjustable Headrest



For Models

- LEON
- Morgan
- Blazer


Hand Bar



For Models


- LEON
- SABER
- SABER TR
- Morgan
- Blazer
- KP-25.2

Controller Cover




Suitable for All Models

USB Charger




Suitable for All Models

Joystick Protector



Suitable for All Models


LED Lighting System A type



For Models

- Blazer


LED Lighting System B type



For Models

- Leon
- Morgan


Lighting System



For Models

- Ergo Stand


Elevating Legrest



For Models

- Leon
- SABER
- Morgan
- Blazer
- KP-25.2


Tie-down Bracket



For Models

- Blazer


Assistant Brake



For Models

- Ergo Stand
- Blazer
- Ergo Nimble
- KP-10.3 CPT


H-shaped Belt



For Models

- Ergo Stand
- SME

Mud Guard



For Models

- Ergo Stand

Specification

Series	Prodcut Name		Seat Width (cm)	Seat Depth (cm)	Front Seat Height (cm)	Backrest Height (cm)	Seat Angle (X°)	Back Angle (X°)	Armrest Height (cm)	Legrest Length (cm)	Power Funtions	Suspension	Front/Middle / Rear Wheel (inch)
Outdoor / Indoor	EVO ALTUS		42-52	40-60	41.5-46.5	51-67	-	-	18-30	40-60	Standing Tilt: 0°-45° Recline: (-5)°-90° Elevating legrest:(-5)°-90° Lift: 0-40cm	4 wheels	14/ - / 9
	MID LECTUS		42-52	40-60	46	51-67	-	-	18-30	40-60	Tilt: 0°-45° Recline: (-5)°-90° Elevating legrest:(-5)°-90° Lift: 0-30cm	6 wheels	8/ 14 / 6
	EVO LECTUS		42-52	40-60	45	51-67	-	-	18-30	40-60	Tilt: 0°-50° Recline: (-5)°-90° Elevating legrest:(-5)°-90° Lift: 0-40cm	4 wheels	14/ - / 9
	LEON	Sling Seat	43/48	41-46	46 (ground to seat pan)	43	2, 5	(-1)-11	18-28	39-51	Tilt: 5°-23° Recline: 0°-32°	4 wheels	9/ -/ 14
		Captain Seat	46	43	51	53	-	0-65	21-31	39-52	-	4 wheels	9/ -/ 14
			51	48		48		0-45					
		Kiss Adult	41-51	41-51	44-49	51-56	0, 3, 6	0-30	20-30	32-46	Tilt: 0°-47° Recline: 0°-80° Elevating legrest: 0°-80° Lift: 0-30cm	4 wheels	9/ -/ 14
		Kiss Junior	31-41	31-41	44-49	41-46			15-25	28-37			
Outdoor	SABER		45	43	53	54	-	5-65	21-27	37-44	-	ISR sys.	9.5/ -/ 16
			51	48									
	SABER TR		51	48	59	54	-	-	21-27	43-52	Tilt: 9°-47° Recline: 0°-25°	ISR sys.	9.5/ -/ 16
Indoor / Outdoor	Ergo Stand		40/44/48	44-50	54	44	-	-	20-28	38-44	Power stand Recline: (-10)°- 30°	-	12 / - / 6
	SME		40/44/48	44-50	52	42	-	-	18-26	36-42	Power stand	-	6 / - / 24
	eFlexx	EFL- 14	38/43/48	40.5/43/45.5	49.5	43	5	(-4)-16	21-29	38-48	-	-	8 / - / 14
		EFL- 20											8 / - / 20
	Morgan	Sling Seat	43/48	41-46	46 (ground to seat pan)	43	2,5	(-1) - 11	18-28	39-51	Tilt: 5°-23° Recline: 0°-32°	6 wheels	8 / 14 / 6
		Captain Seat	46	43	51	53	-	0 - 65	21-31	39-52	-	6 wheels	8 / 14 / 8
			51	48		48		0 - 45					
		Kiss Adult	41-51	41-51	44-49	51-56	0,3,6	0 - 30	20-30	32-46	Tilt: 0°-47° Recline: 0°-80° Elevating legrest: 0°-80° Lift: 0-30cm	6 wheels	8 / 14 / 6
		Kiss Junior	31-41	31-41	44-49	41-46			15-25	28-37			
	Blazer	Sling Seat	43/48	41-46	44 (ground to seat pan)	43	-	(-1) - 14	18-28	32-46	Tilt:4°-25°	rear	9 / - / 14
		Captain Seat	46	43	54	53	-	0 - 65	21-31	37-48	-	rear	9 / - / 14
			51	48									
	KP-25.2		36	43	50	43	-	-	23-28.5	39-51	-	-	8 / - / 12
			41	43									
			46/51	46									
Ergo Nimble		44	43	54	48	-	-	21-27	38-46	-	-	8/ - / 12	
KP-10.3 CPT		45	46	52	44	-	15	23-27	38-46	-	-	8/ - / 12	
SP-100		41/46	42	49	38	6	6	22	37-43	-	-	8/ - / 16	

Con-troller	Motor	Battery	Min. Turning Radius (cm)	Min. Ground Clearance (cm)	Max. Safe Slope (X°)	Max. Obstacle Height (cm)	Max. Speed (km/h)	Max. Driving Range (km)	Max. User Weight (kg)	w/ Batteries / w/o Batteries	Overall L x W x H (cm)
Rnet	320W	50Ah*2/ 80Ah*2	75	7.5	10	7	6,10,12	50Ah=25/ 80Ah=45	140	215(80Ah) / 165	120*61*112
Rnet	350W	50Ah*2/ 80Ah*2	62	9	10	7.5	6,10,12	50Ah=24/ 80Ah=36	140	190(80Ah) / 133	110*61.5*104
Rnet	320W	80Ah*2	75	7	10	7	6,10	40	140	210 / 160	120*62*112
VR2	450W	50Ah*2/ 80Ah*2	78	6.5	8	5	6, 10, 12	30 (50Ah)/ 45 (80Ah)	140	108/ 78	106*62*130
VR2	450W	50Ah*2/ 80Ah*2	78	6.5	8	5	6, 10, 12	30 (50Ah)/ 45 (80Ah)	140	100/ 70	106*62*123
VR2/ Rnet	450W	50Ah*2/ 80Ah*2	90	6.5	6	5	10	25 (50Ah)	140	Min. 109/ Min. 79	117*65*124
			85					25 (50Ah)	100	Min. 106/ Min. 76	105*65*110
VR2	400W	50Ah*2/ 80Ah*2	88	10	10	8	13	38 (50Ah)/ 57 (80Ah)	135	102/ 73	117.5*65*120
VR2	400W	50Ah*2	97	10	10	8	13	38	135	113/ 82	120*65*126
VR2	300W	36Ah*2/ 50Ah*2	67	6.5	10	5	9	28 (50Ah)	135	89/57	102*63*98
-	-	2.9Ah*2	-	-	-	-	-	-	120	28.6/25.8	104*(59-67)*(96-143)
Linx	250W	26.1/11.5	90	7	6	2.5	8	20	100	33(STD) /29 35.5 (STD) /31.5	110*65*85 117.5*66*85
VR2	450W	50Ah*2	64	8	8	5	6, 10	25	140	116/86	115*60*130
VR2	450W	50Ah*2	64	8	8	5	6,10	25	140	108/78	115*60*123
VR2	450W	50Ah*2	72.5	8	6	5	10	25	140	117/87	122*61*103
			67.5					25	100	114/84	117*61*103
VR2	300W	50Ah*2	82	7	8	6	10	30	140	94/64	107*60*93
VR2	300W	50Ah*2	82	7	8	6	10	30	140	100/70	107*63*128
VR2	200W	22Ah*2 36Ah*2/ 50Ah*2 36Ah*2/ 50Ah*2	82	6	10	5	6	38(50Ah)	115 (135 pne.)	70/35	110*(54-69)*96
VR2	200W	22Ah*2	75	9	8	3	6	12	120	50/35	103*62*94
VR2	200W	22Ah*2	75	9	8	3	6	12	120	Min.53.3/ Min.37.5	103*61*92
Linx	100W	22Ah*2	81	7	6	4	6	15	100	43/29	99*(56/61)*86



A Better Wheelchair, A Better Fit

Karma Healthcare Ltd.

9th floor, Bhakta Towers, KB 22,
Salt Lake, Kolkata - 700 098

Phone: +91-33-2335 7062

NORTH

DELHI : 9350834010 / 9821113831,
JAIPUR : 8890050487,
PUNJAB : 9821113831,
LUCKNOW : 7439020124,
HARYANA : 9891380054,
RAJASTHAN : 8890050487

EAST

PATNA : 9386253881,
BHUBANESWAR : 9438335512,
GUWAHATI : 7439020144,
KOLKATA : 033-23357062

WEST

PUNE & SOUTH MAHARASTRA : 8983501408,
AHMEDABAD : 022 2854 3924,
MUMBAI : 022 2854 3924,
NAGPUR : 9146303373,
INDORE : 022 2854 3924,
GOA : 9769050727

SOUTH

CHENNAI : 9381794124 / 9362230525,
VIJAYAWADA : 9397602122,
BANGALORE : 9788408666,
COCHIN : 9944094587,
COIMBATORE : 9362230525,
HYDERABAD : 9397602122

Email: service@karmahealthcare.com

<http://www.karmahealthcare.com/>



ISO 9001